



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Salbutamol
---------------	------------

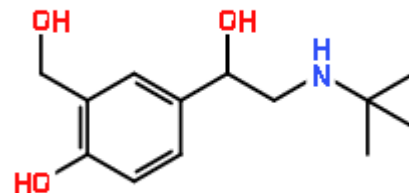

CHROMATOGRAM

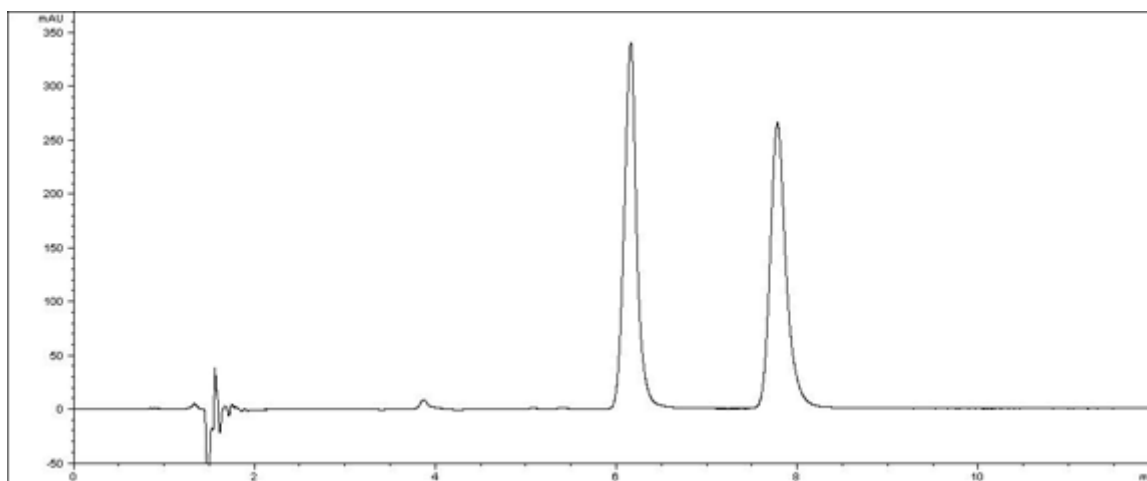

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	12.5mM 2-amino-1-butanol in methanol / tetrahydrofuran / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.15	2.34	1 - 2	6.04
2	7.78	3.22		