



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Zopiclone
---------------	-----------

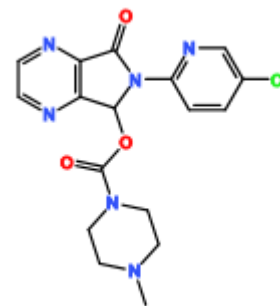

CHROMATOGRAM

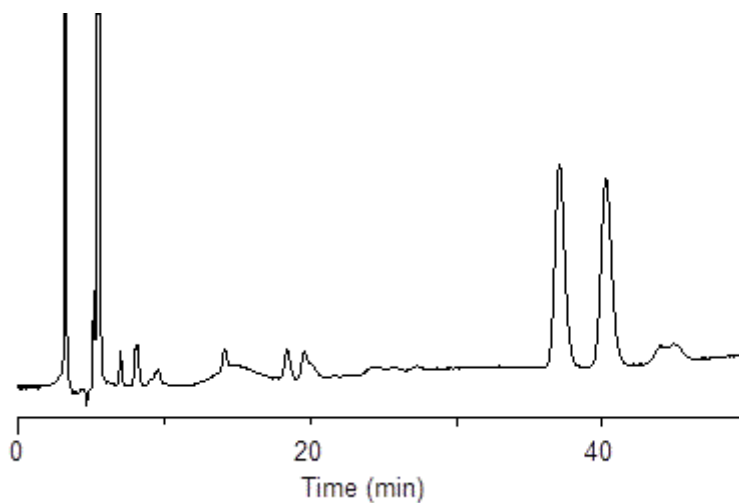

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IE 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 80 / 20 / 0.1
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	37.1	11.4	1 - 2	2.49
2	40.3	12.4		