


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Terbutaline

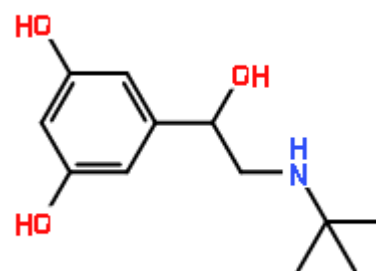

## CHROMATOGRAM

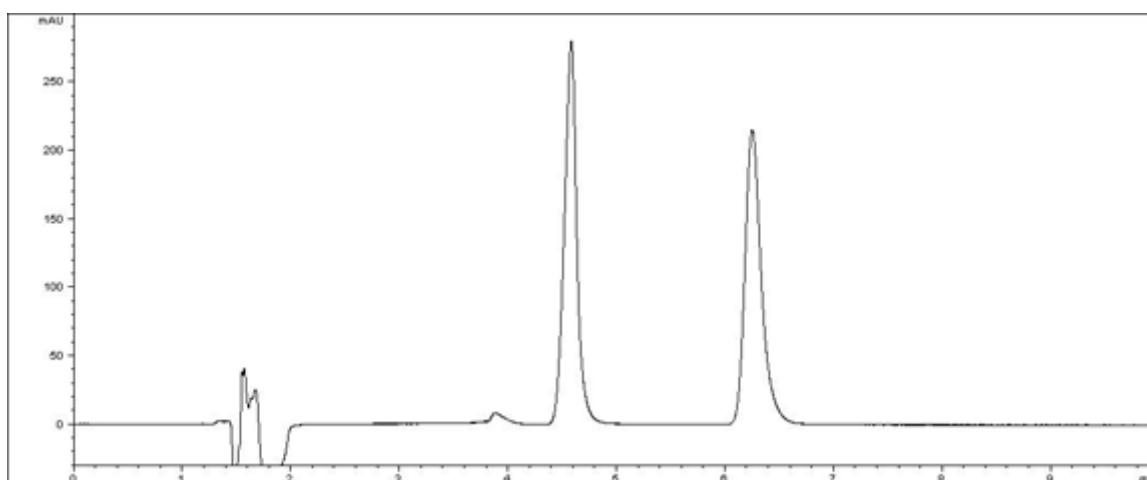

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	12.5mM formic acid + 25mM diethylamine in methanol / tetrahydrofuran / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.58	1.69	1 - 2	7.24
2	6.25	2.67		