



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	trans-4-Hydroxy-DL-proline
---------------	----------------------------

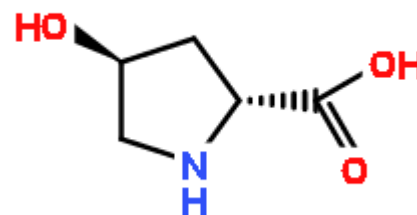

CHROMATOGRAM

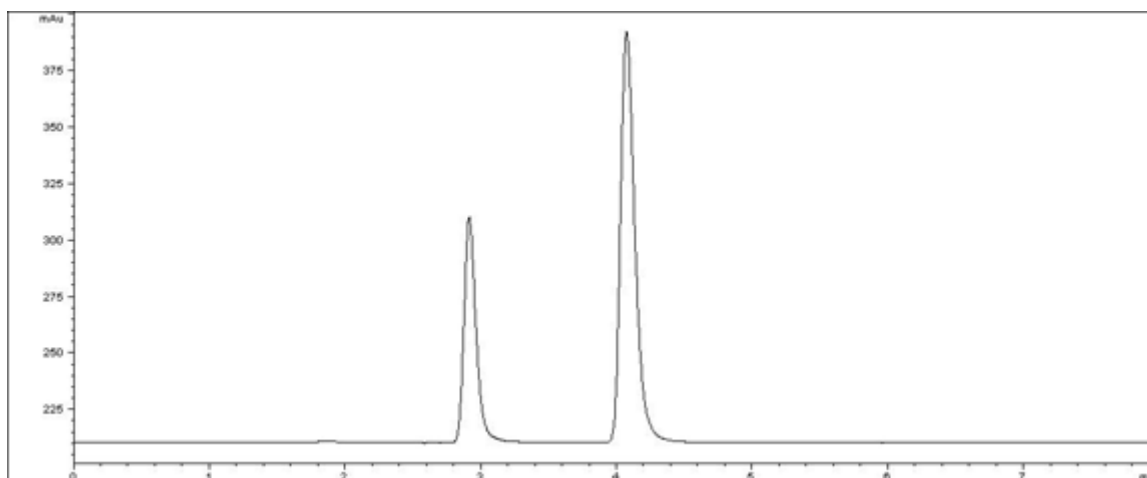

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(-) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / acetonitrile / water = 49 / 49 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0015
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	D/L

No.	Rt	K'	Alpha Value	Resolution
1	2.91	0.64	1 - 2	6.81
2	4.07	1.29		