


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ethiazide

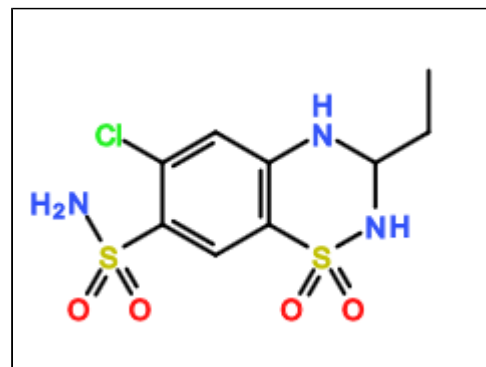

## CHROMATOGRAM

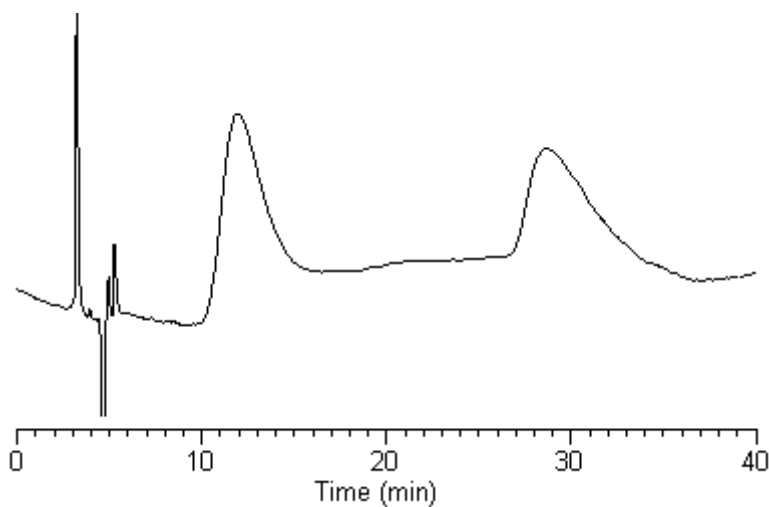

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IE 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 80 / 20 / 0.1
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	12	3	2.83	3.41
2	28.6	8.5		