



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	(±)-2,3-O-Benzylidene threitol
---------------	--------------------------------

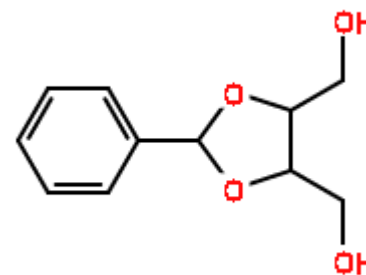

## CHROMATOGRAM

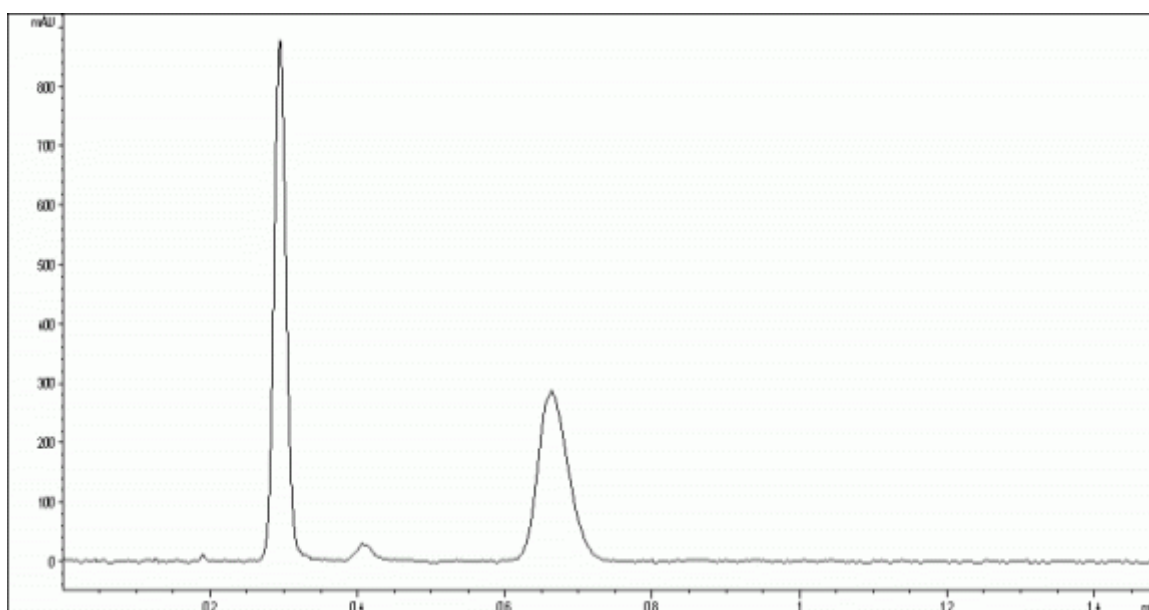

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA-3 50 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / ethanol = 80 / 20
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	5	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.3	1.57	1 - 2	6.9
2	0.66	4.77		