



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	(4R,5S/4S,5R)-1,5-Dimethyl-4-phenyl-2-imidazolidinone
---------------	---

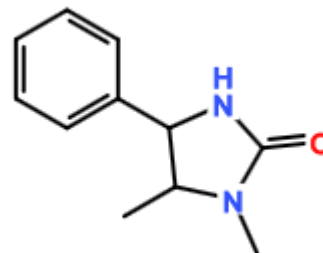

CHROMATOGRAM

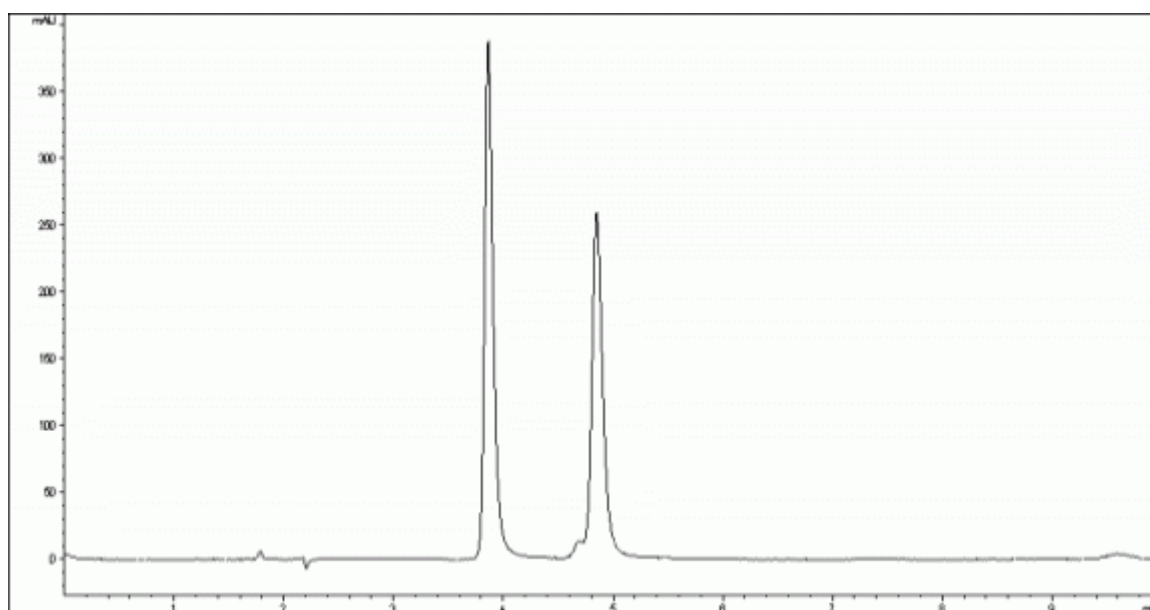

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol = 80 / 20
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.0025
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.86	1.27	1 - 2	6.43
2	4.85	1.85	1.46	