



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	5-Benzylhydantoin
---------------	-------------------

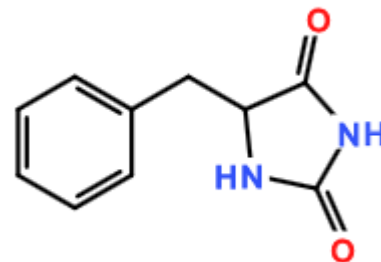

## CHROMATOGRAM

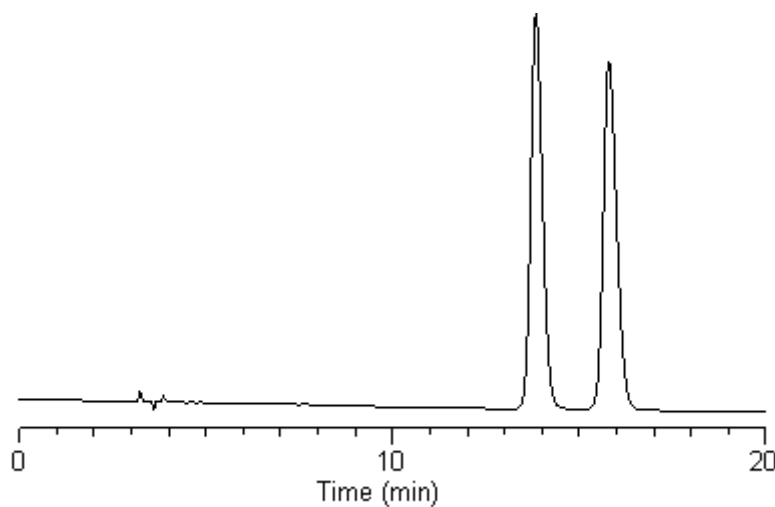

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IF 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol = 70 / 30
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 220 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	13.8	3.6	1 - 2	3.19
2	15.8	4.3		