


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Mitiglinide

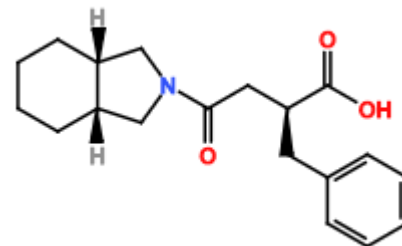

CHROMATOGRAM

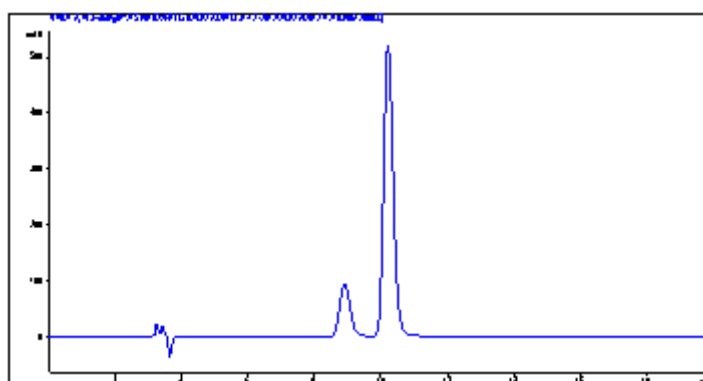

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methyl tert-butyl ether / acetonitrile / trifluoroacetic acid = 80 / 20 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 215 nm		
Elution Order	Enantiomer, API		

No.	Rt	K'	Alpha Value	Resolution
1	8.9	1.97	1 - 2	2.38
2	10.2	2.40		