



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Ramelteon (racemic)
---------------	---------------------

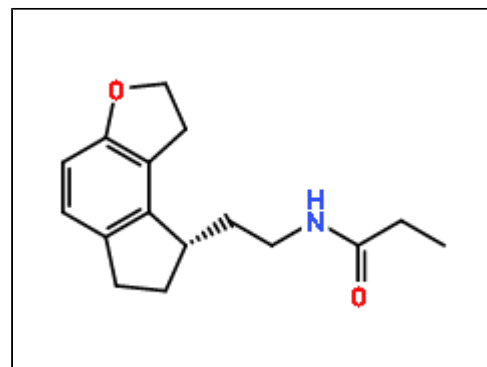

CHROMATOGRAM

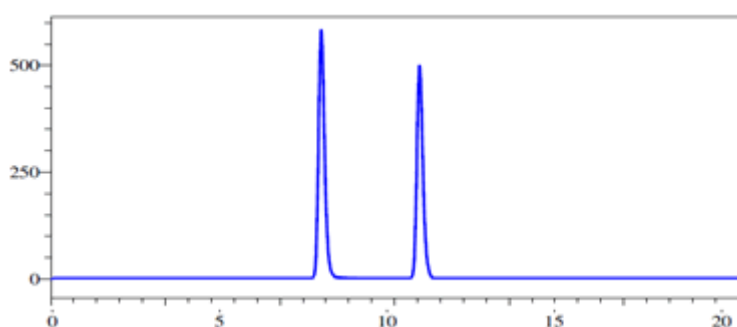

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ID 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 80 / 20 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.02
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 285 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	8.12	1.71	1 - 2	6.14
2	10.94	2.65		