


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Rasagiline

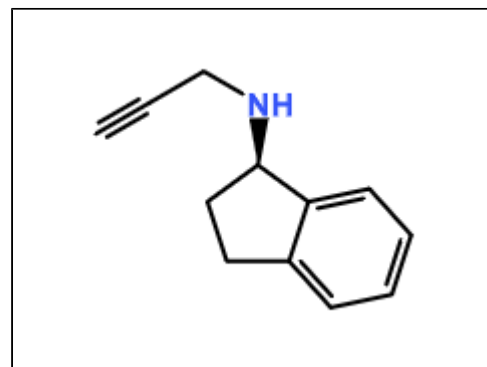

CHROMATOGRAM

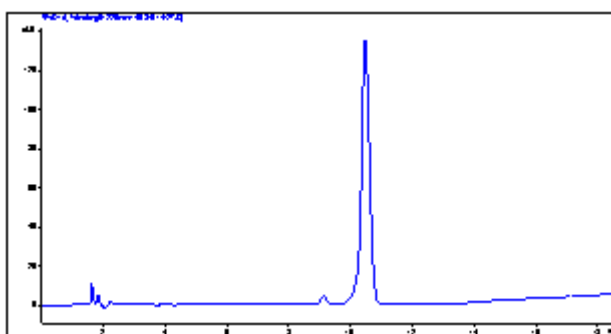

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IF-3 150 x 4.6mm / 3µm
Mobile Phase Composition	10mM ammonium bicarbonate aq. / methanol / acetonitrile = 50 / 40 / 10
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 220 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.14	4.08	1 - 2	1.18
2	10.48	4.82		3.21