


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Rebamipide

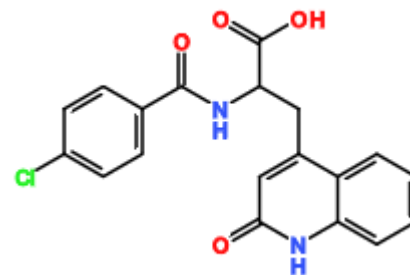

CHROMATOGRAM

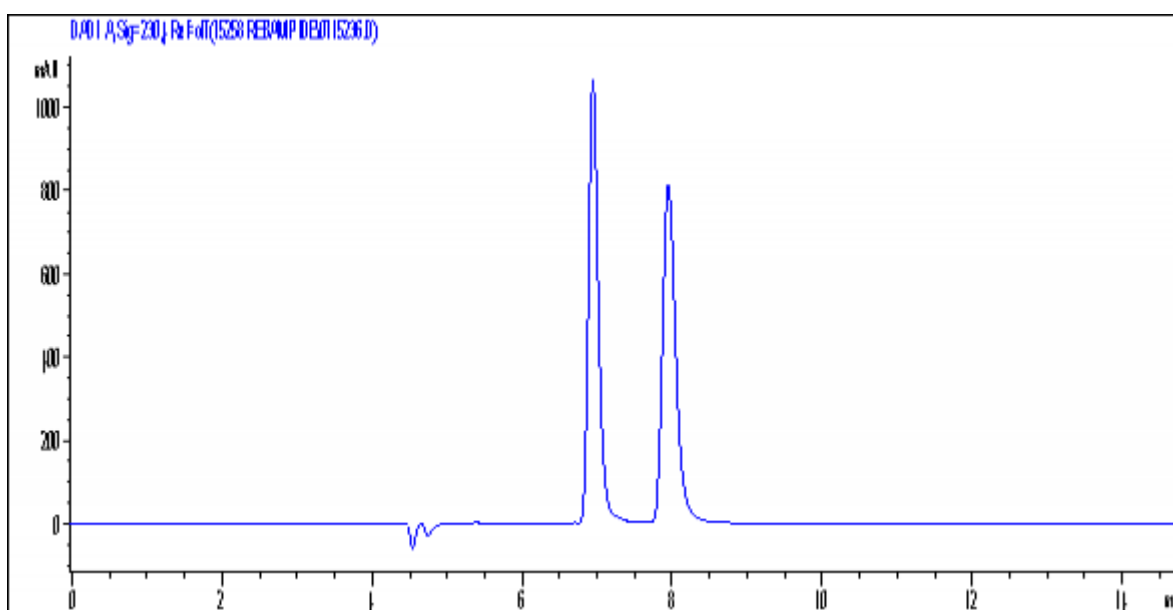

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IE-3 250 x 4.6mm / 3µm
Mobile Phase Composition	acetonitrile / methanol / trifluoroacetic acid = 85 / 15 / 0.3
Chromatographic Mode	Normal Phase Secondary screening solvent (Immobilised CSP)

Flow rate (ml/min)	0.7	Injection amount (mg)	0.002
Temperature (°C)	40	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.94	0.62	1 - 2	3.74
2	7.95	0.86		