



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Valganciclovir
---------------	----------------

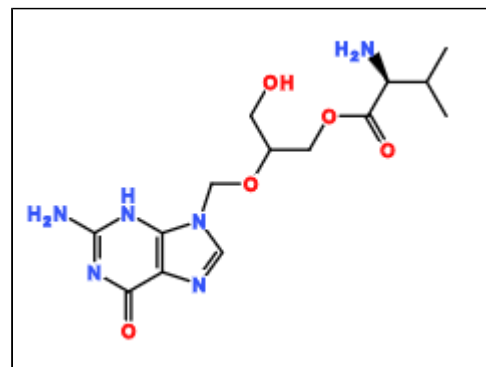

CHROMATOGRAM

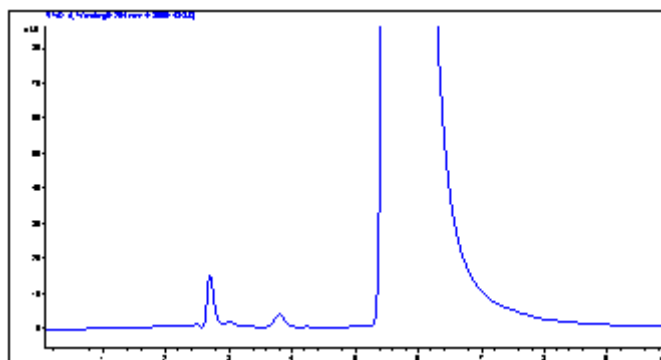

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CROWNPAK® CR-I(+) 150 x 3mm / 5µm
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 95 / 5
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.4	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	D, L

No.	Rt	K'	Alpha Value	Resolution
1	3.8	0.77	1 - 2	2.46
2	6.2	1.88		4