



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Ethotoin
---------------	----------

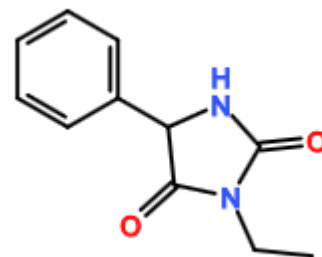

CHROMATOGRAM

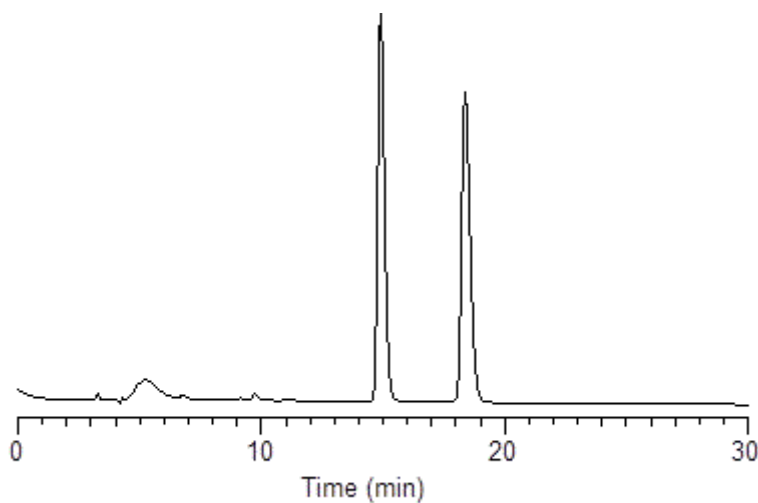

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol = 90 / 10
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 220 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	14.9	4	1 - 2	6.07
2	18.4	5.1		