



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Garenoxacin mesylate
---------------	----------------------

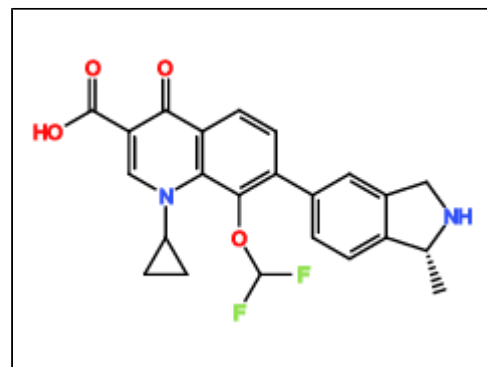

## CHROMATOGRAM

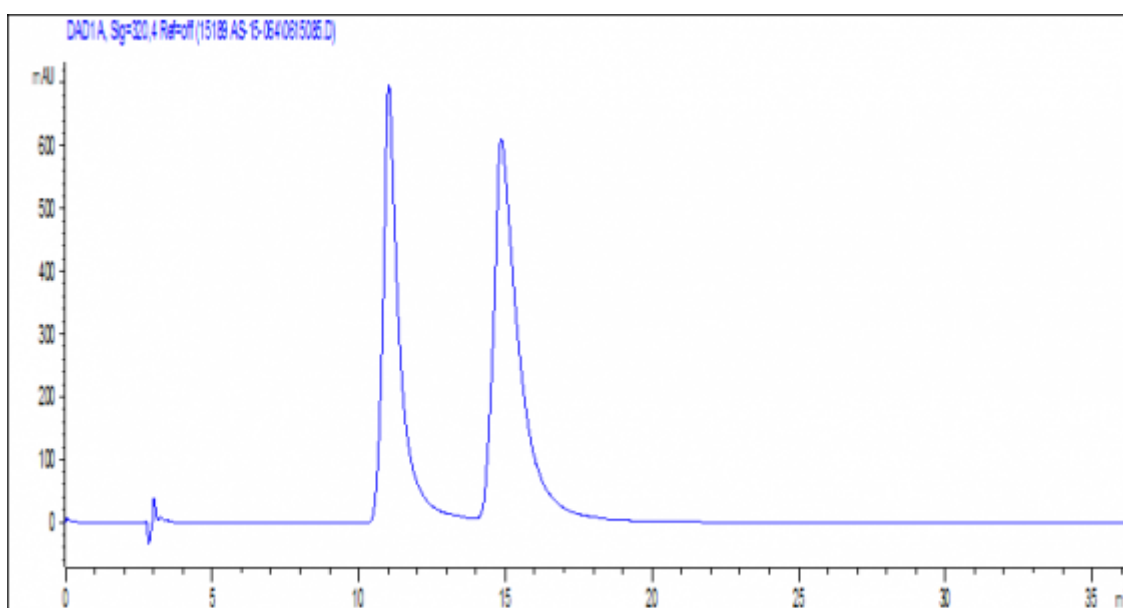

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	100mM formic acid + 25mM diethylamine in methanol / tetrahydrofuran / water / acetic acid = 49 / 49 / 2 / 0.5
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.02
Temperature (°C)	15	Pressure (MPa)	
Detection	UV 320 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	11.02	2.06	1 - 2	3.5
2	14.85	3.13		