



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Ramelteon (racemic)
---------------	---------------------

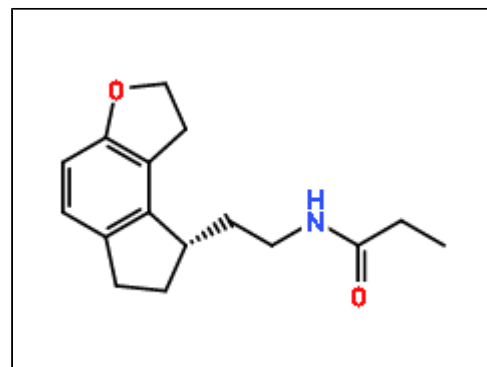

CHROMATOGRAM

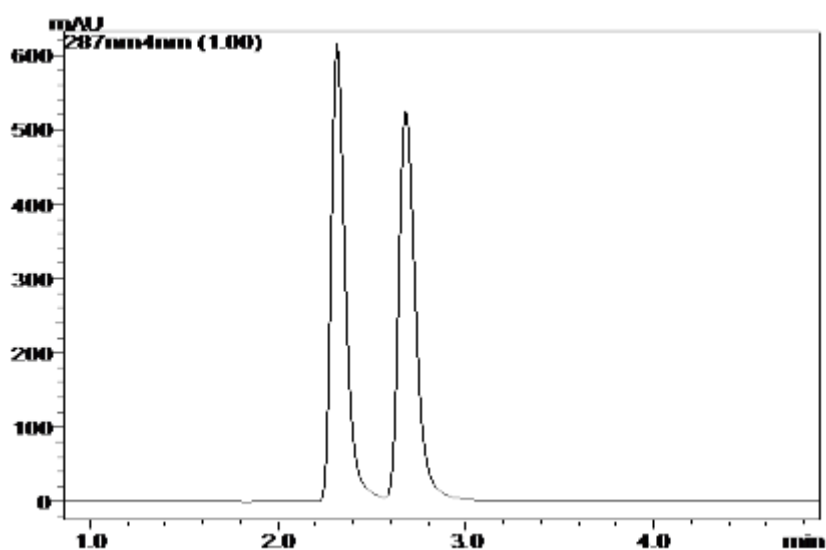

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AD-3R 150 x 4.6mm / 3µm
Mobile Phase Composition	100mM ammonium bicarbonate aq. / methanol = 5 / 95
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 287 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.31	0.28	1 - 2	2.21
2	2.68	0.49		