


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Zopiclone

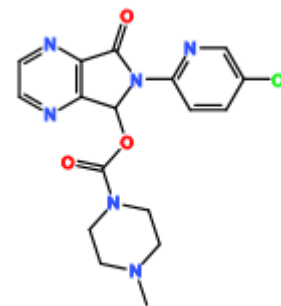

CHROMATOGRAM

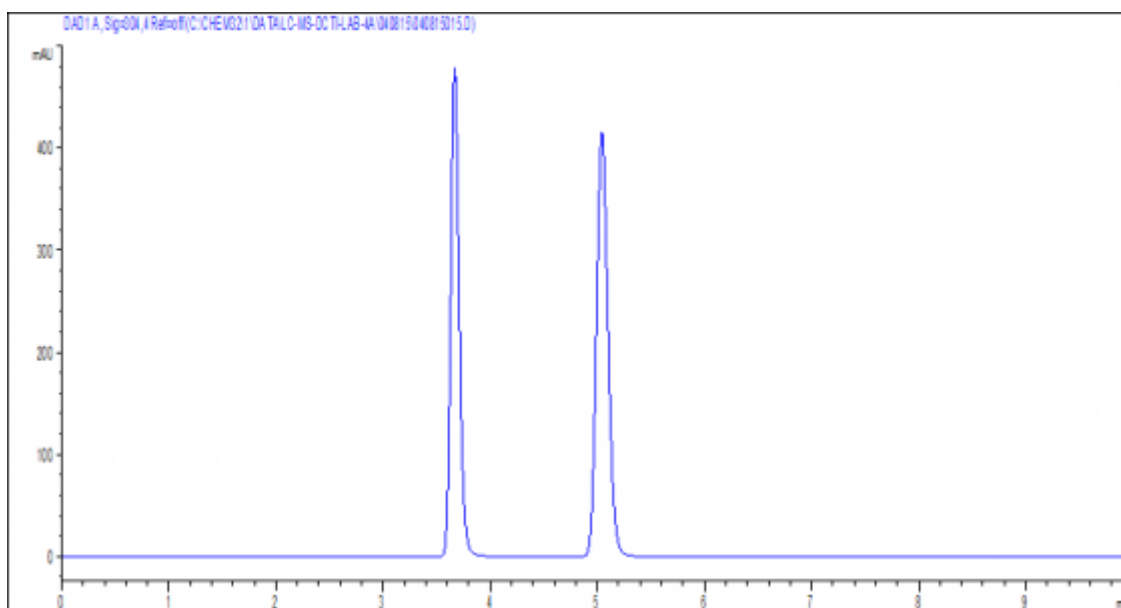

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC-3 150 x 4.6mm / 3µm
Mobile Phase Composition	100mM ammonium bicarbonate aq. / acetonitrile = 5 / 95
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	40	Pressure (MPa)	
Detection	UV 304 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.67	1.04	1 - 2	1.73
2	5.04	1.8		8.07