



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Ranolazine
---------------	------------

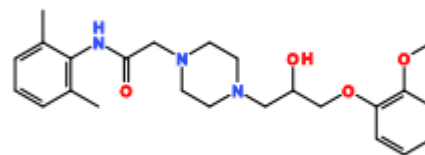

CHROMATOGRAM

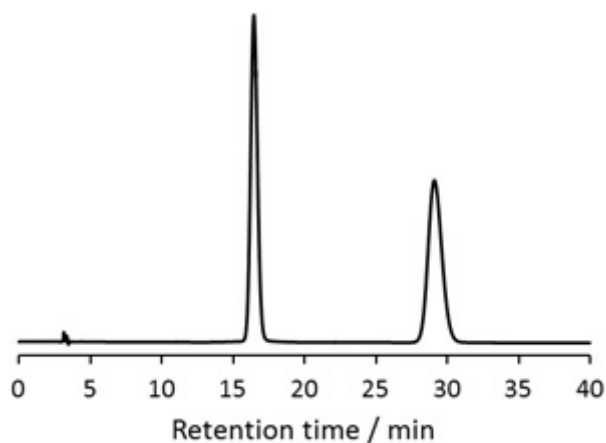

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IG 250 x 4.6mm / 5µm
Mobile Phase Composition	hexane / ethanol / acetic acid / diethylamine = 50 / 50 / 0.1 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	230	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	16.46	4.49	1 - 2	1.94
2	29.1	8.7		9.34