


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cetirizine

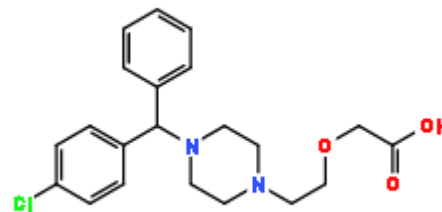

## CHROMATOGRAM

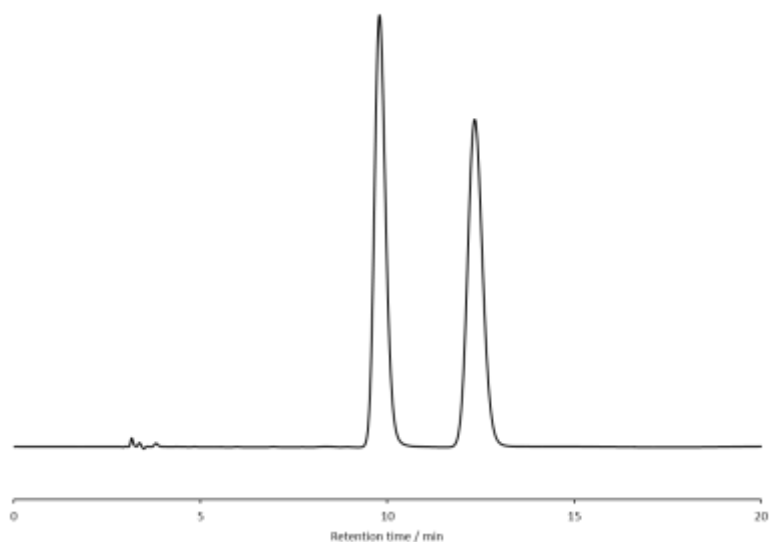

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IG 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol / acetic acid / diethylamine = 50 / 50 / 0.1 / 0.1
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.79	2.26	1 - 2 1.38	3.87
2	12.34	3.11		