



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Brinzolamide
---------------	--------------

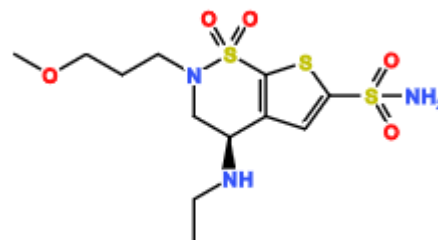

CHROMATOGRAM

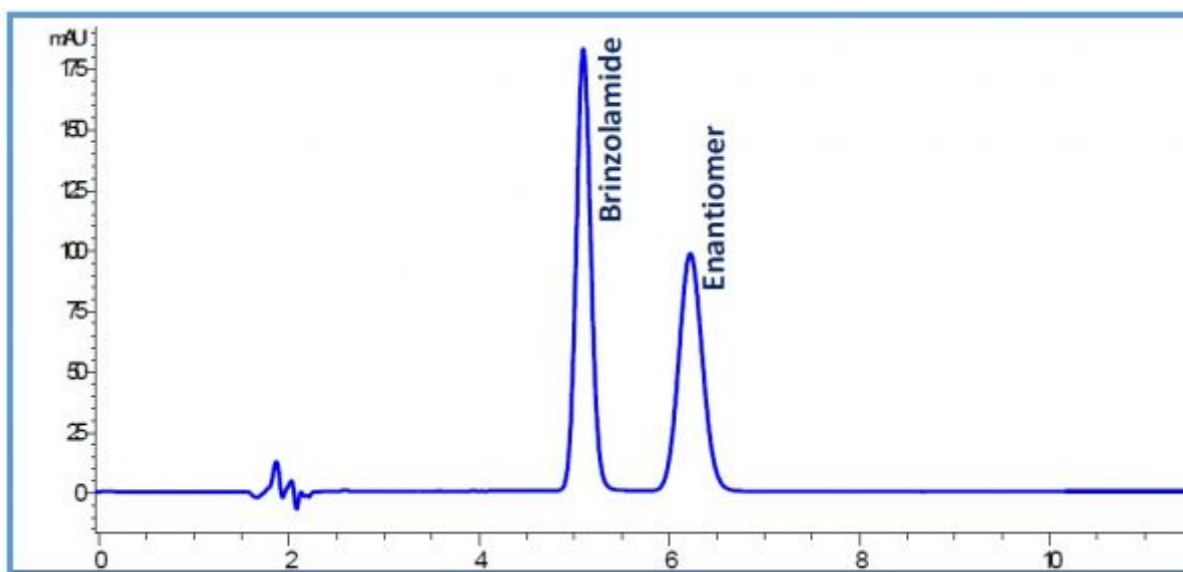

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IG-3 150 x 4.6mm / 3µm
Mobile Phase Composition	10mM ammonium bicarbonate aq. / acetonitrile = 70 / 30
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.1	1.73	1 - 2	2.99
2	6.2	2.32		