



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	4-Phenyl-2-oxazolidinone
---------------	--------------------------

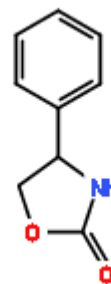

CHROMATOGRAM

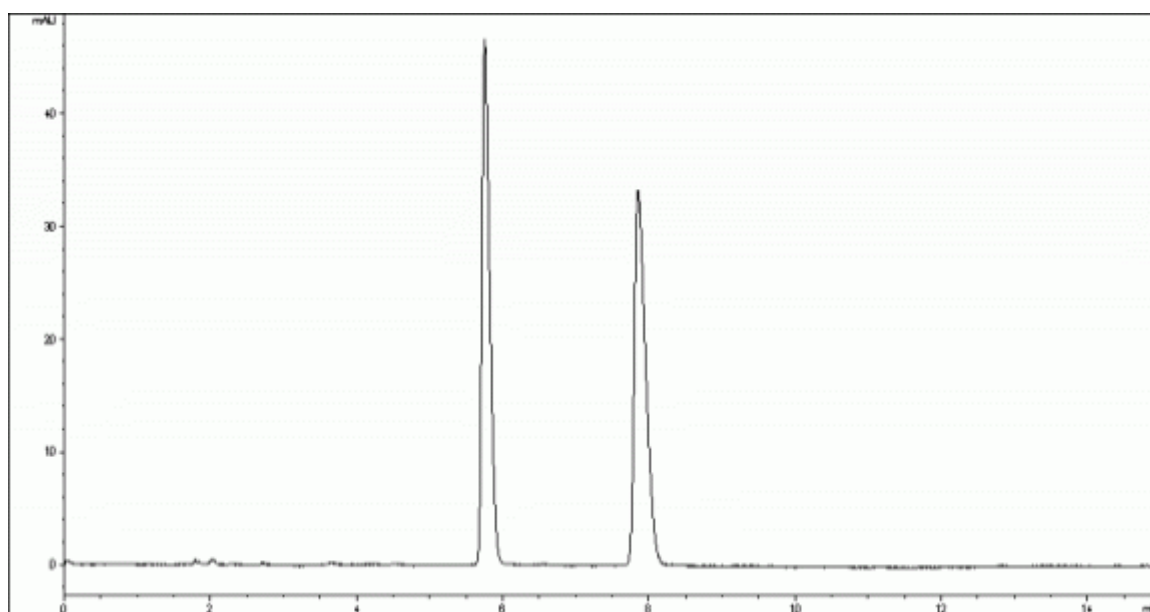

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC-3 150 x 4.6mm / 3µm
Mobile Phase Composition	dichloromethane / ethanol = 98 / 2
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.02
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.75	2.36	1 - 2	8.79
2	7.85	3.59		