


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

DL-Homophenylalanine

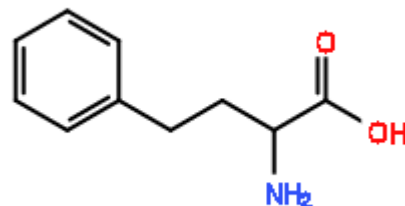

CHROMATOGRAM

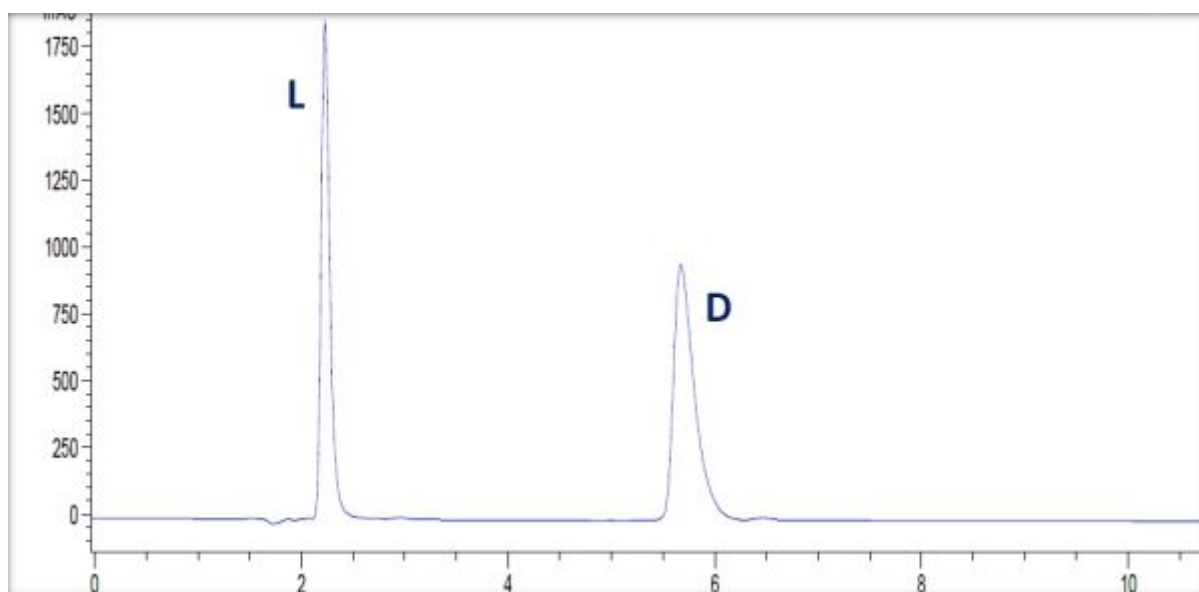

ANALYTICAL REPORT

Column Trade Name	CROWNPAK® CR-I(-)		
Column Size (Particle Size)	3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.0 / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.4	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 215 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.2		1 - 2	13.80
2	5.7			