



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	α -Amino- ϵ -caprolactam
---------------	--

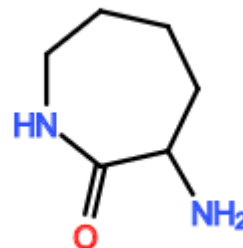

CHROMATOGRAM

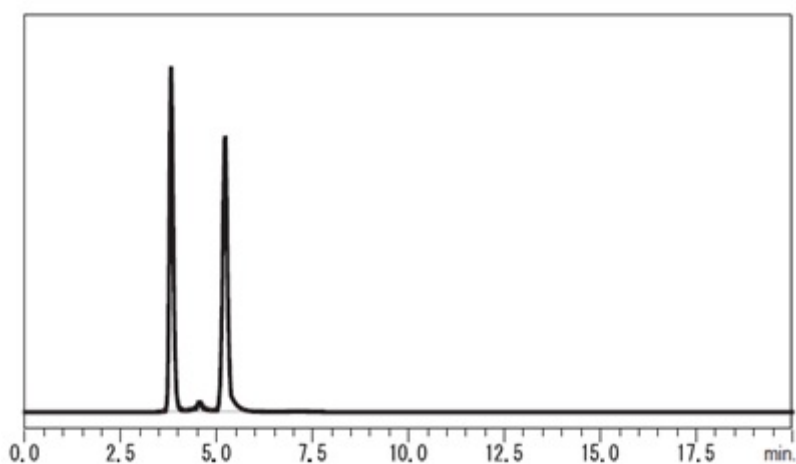

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CROWNPAK® CR-I(+) 150 x 3mm / 5 μ m
Mobile Phase Composition	perchloric acid aq. pH1.5 / acetonitrile = 80 / 20
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.2	Injection amount (mg)	0.001
Temperature (°C)	25	Pressure (MPa)	
Detection	200	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1			1 - 2	