



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Serine
---------------	--------

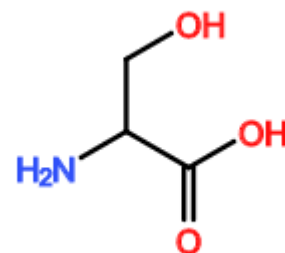

CHROMATOGRAM

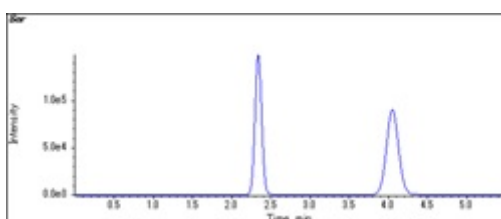

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CROWNPAK® CR-I(+) 150 x 3mm / 5µm
Mobile Phase Composition	acetonitrile / water / trifluoroacetic acid = 96 / 4 / 0.5
Chromatographic Mode	HILIC

Flow rate (ml/min)	0.4	Injection amount (mg)	
Temperature (°C)	30	Pressure (MPa)	
Detection	TOFMS	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	2.34		1 - 2	7.76
2	4.05			