


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Verapamil

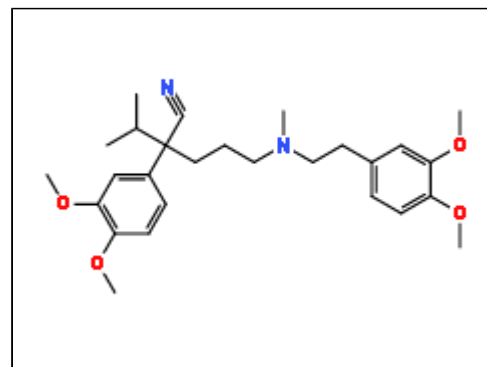

## CHROMATOGRAM

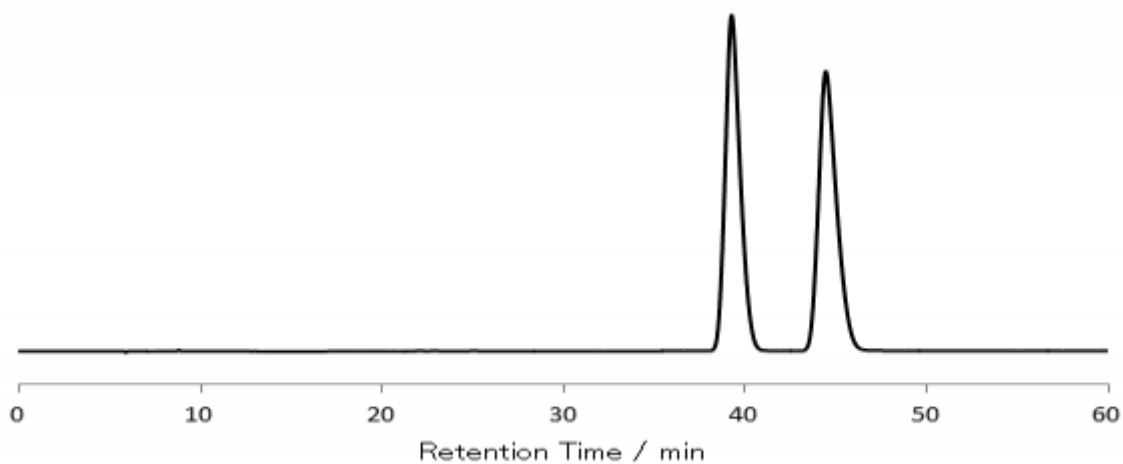

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB N-5 250 x 4.6mm / 5µm
Mobile Phase Composition	0.1M Potassium hexafluorophosphate aq. / acetonitrile = 60 / 40
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	39.2	5.6	1 - 2	1.16
2	44.4	6.5		