


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Clorprenaline

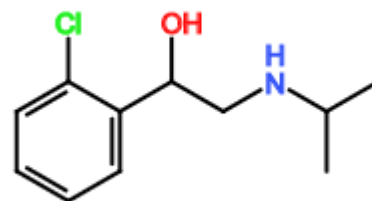

CHROMATOGRAM

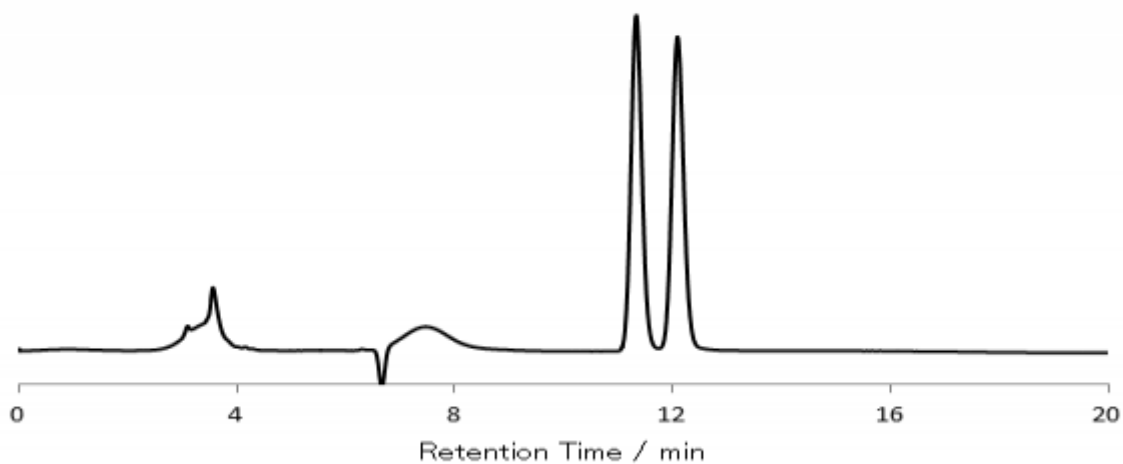

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IG 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 98 / 2 / 0.1
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	11.3	2.8	1 - 2 1.09	1.93
2	12.1	3		