



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	gamma-Phenyl-gamma-butyrolactone
---------------	----------------------------------

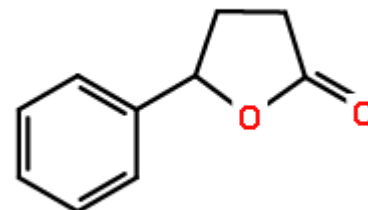

CHROMATOGRAM

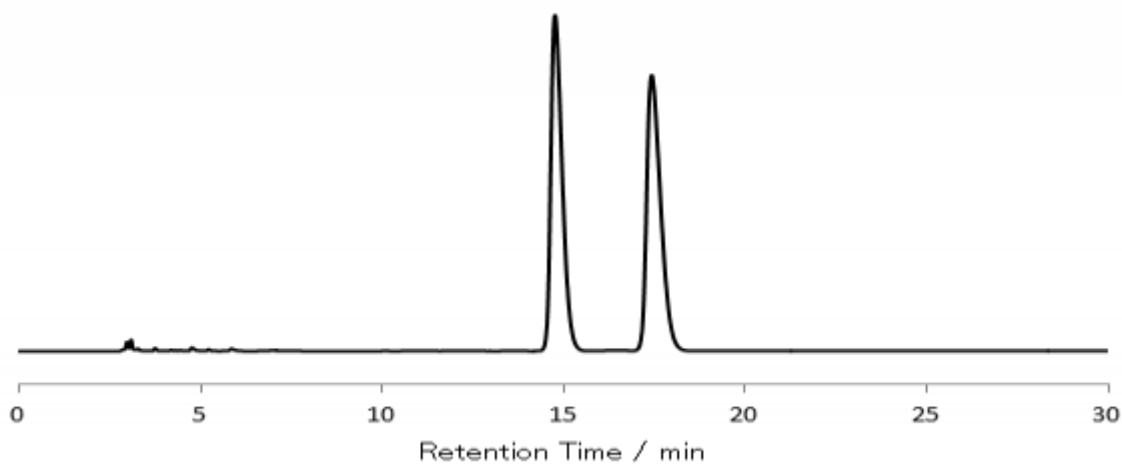

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IH 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol = 80 / 20
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	14.8	3.9	1 - 2	4.19
2	17.4	4.8		