



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Flurbiprofen
---------------	--------------

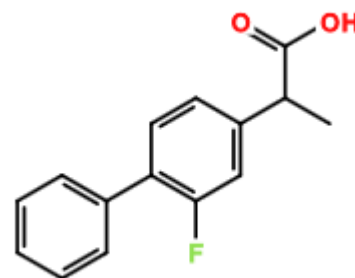

CHROMATOGRAM

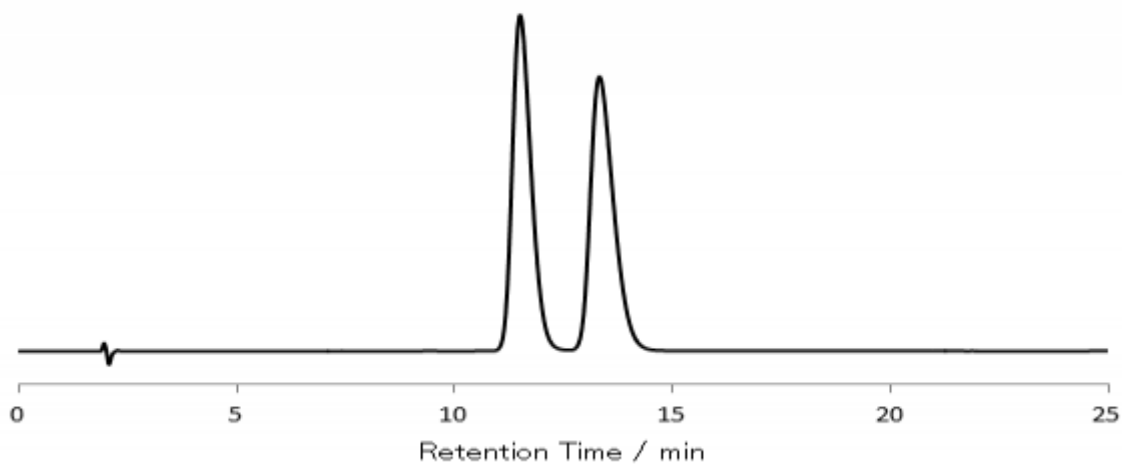

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IJ 150 x 4.6mm / 5µm
Mobile Phase Composition	formic acid aq. pH2.0 / methanol = 20 / 80
Chromatographic Mode	Reverse Phase

Flow rate (ml/min)	1.0	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 245 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	11.5	4.75	1.19	2.07
2	13.3	5.66		