



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Mandelic acid methyl ester
---------------	----------------------------

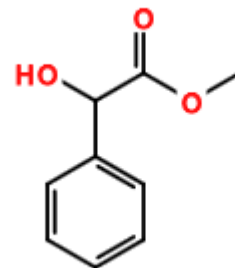

## CHROMATOGRAM

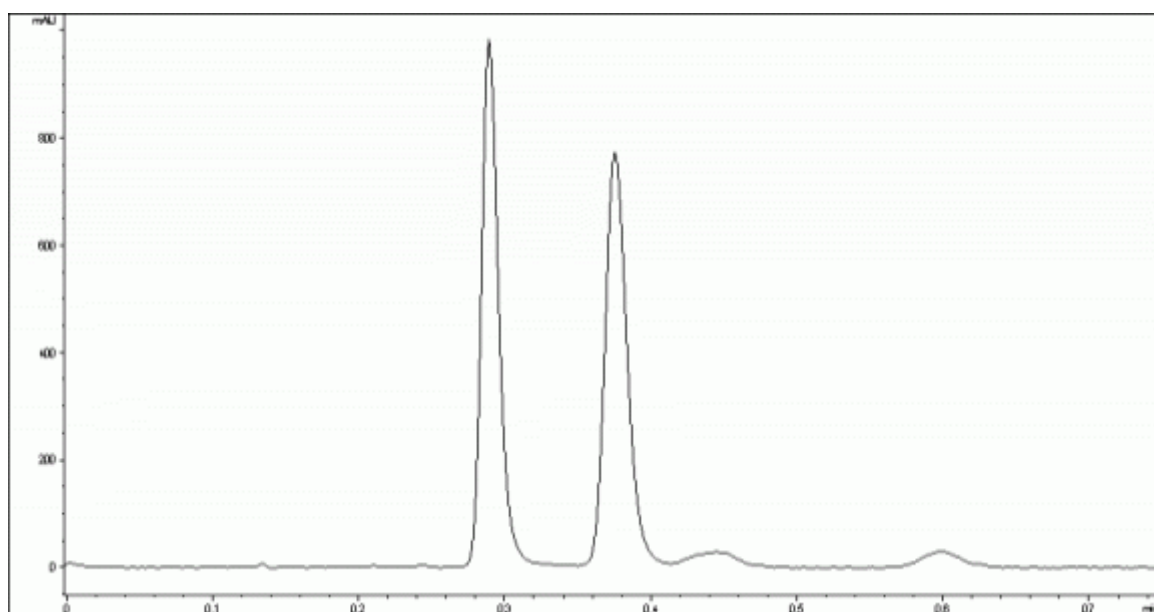

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB-3 50 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10
Chromatographic Mode	Normal Phase

Flow rate (ml/min)	5	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 210 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	0.29	1.33	1 - 2	3.58
2	0.38	2.02		