



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-benzoyl-DL-Leucine
---------------	----------------------

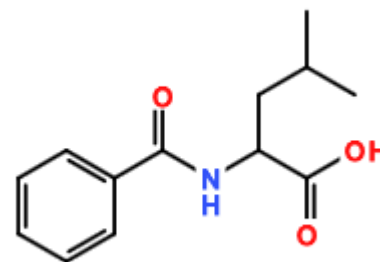

CHROMATOGRAM

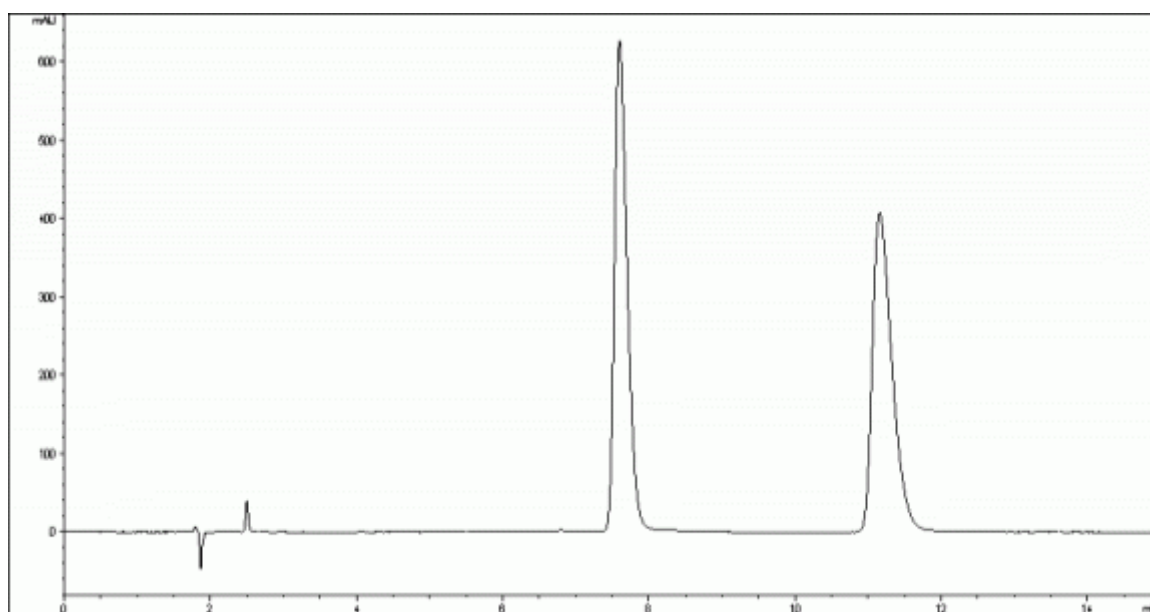

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / tetrahydrofuran / trifluoroacetic acid = 85 / 15 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.6	3.44	1 - 2	1.61
2	11.16	5.53		8.9