



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-FMOC-DL-Leucine
---------------	-------------------

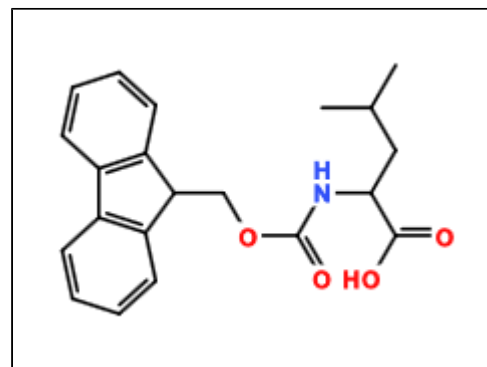

CHROMATOGRAM

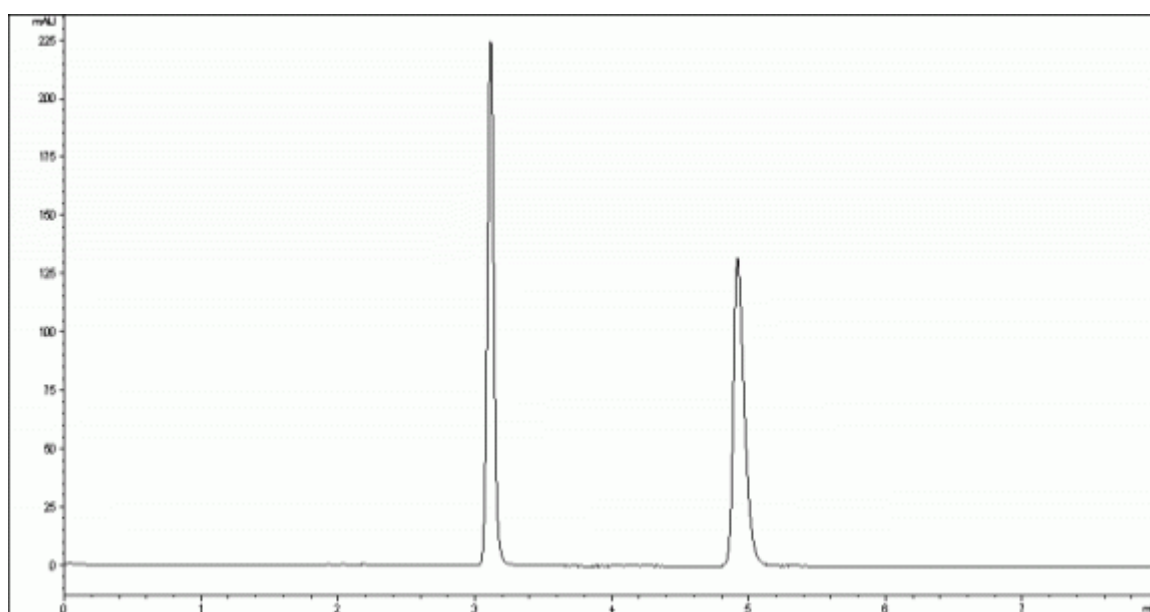

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / chloroform / ethanol / trifluoroacetic acid = 68 / 30 / 2 / 0.1
Chromatographic Mode	Normal Phase Secondary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.12	0.72	1 - 2	16.06
2	4.92	1.71		