



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Tropic acid
---------------	-------------

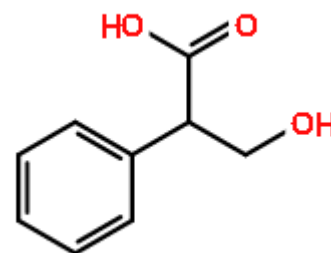

CHROMATOGRAM

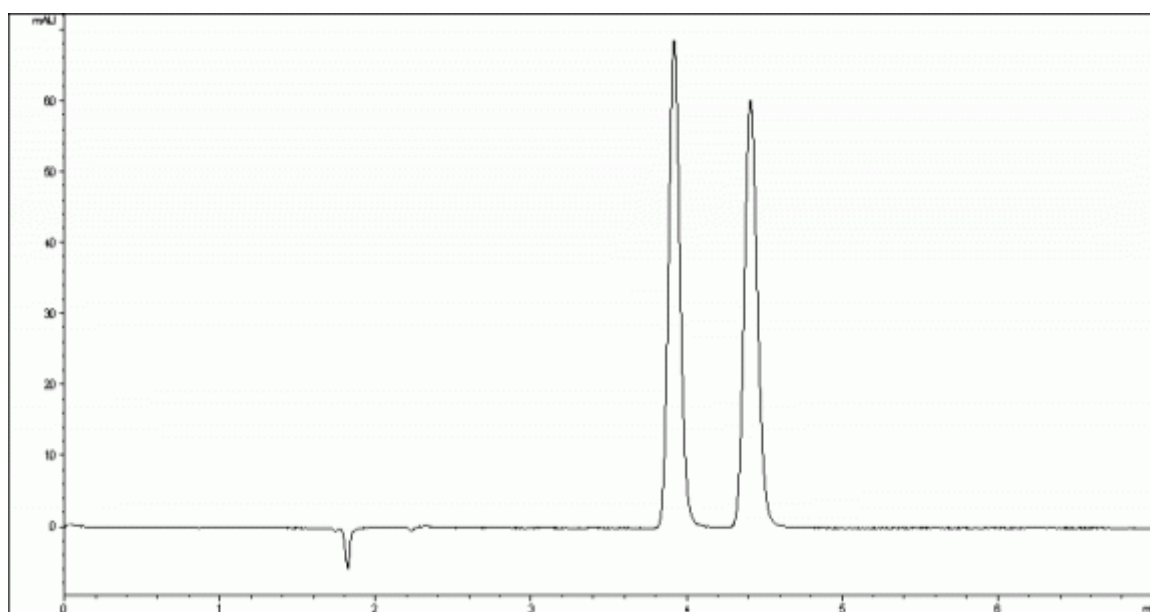

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC-3 150 x 4.6mm / 3µm
Mobile Phase Composition	dichloromethane / ethanol / trifluoroacetic acid = 98 / 2 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.02
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	3.92	1.29	1 - 2	3.44
2	4.41	1.58		