


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Warfarin

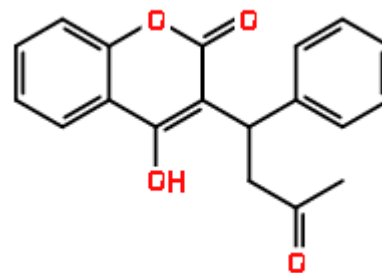

CHROMATOGRAM

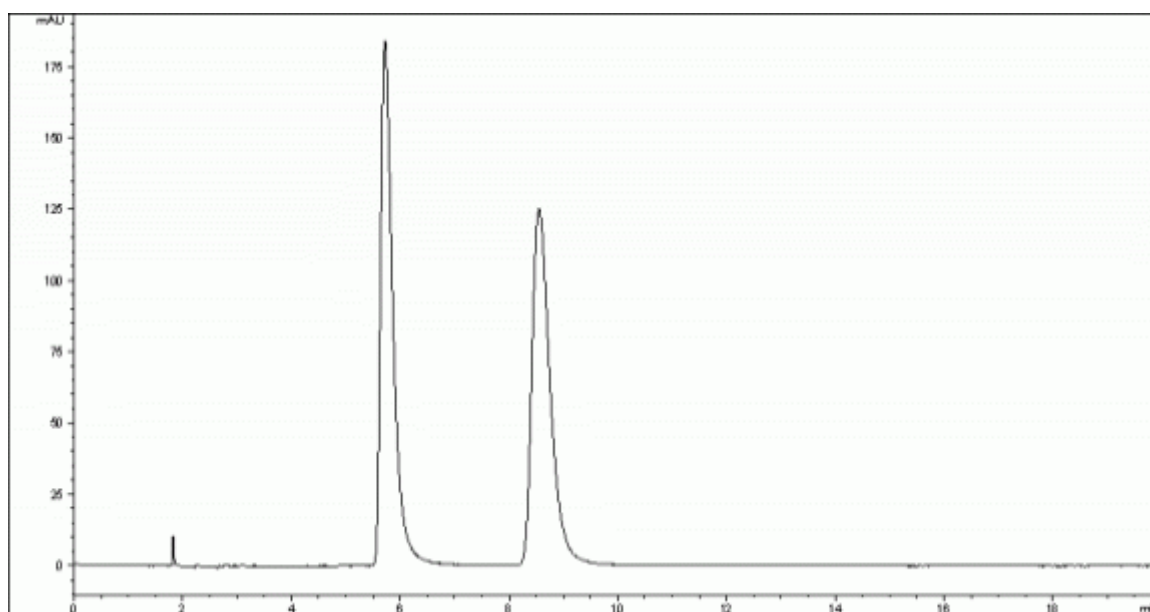

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / ethyl acetate / trifluoroacetic acid = 70 / 30 / 0.1
Chromatographic Mode	Normal Phase Secondary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.72	2.15	1 - 2	5.79
2	8.55	3.71		