



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	4-Benzoyloxy-2-azetidinone
---------------	----------------------------

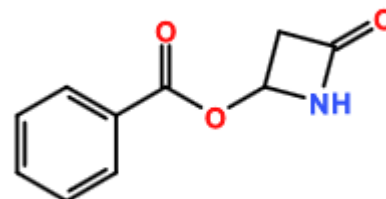

CHROMATOGRAM

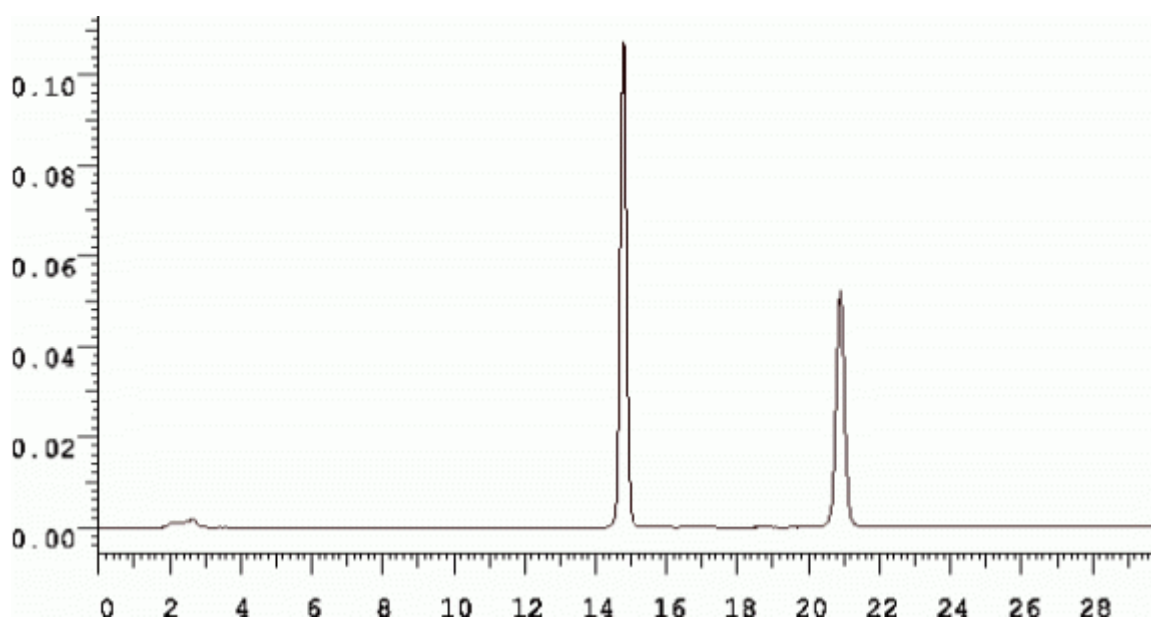

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / dichloromethane = 25 / 75
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	14.8	3.93	1 - 2	16.45
2	20.89	5.96		