



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	beta-Alanyl-DL-leucine
---------------	------------------------

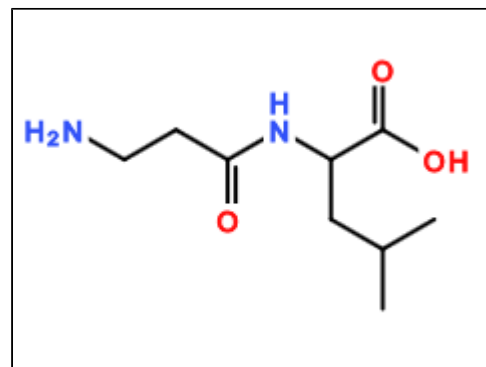

CHROMATOGRAM

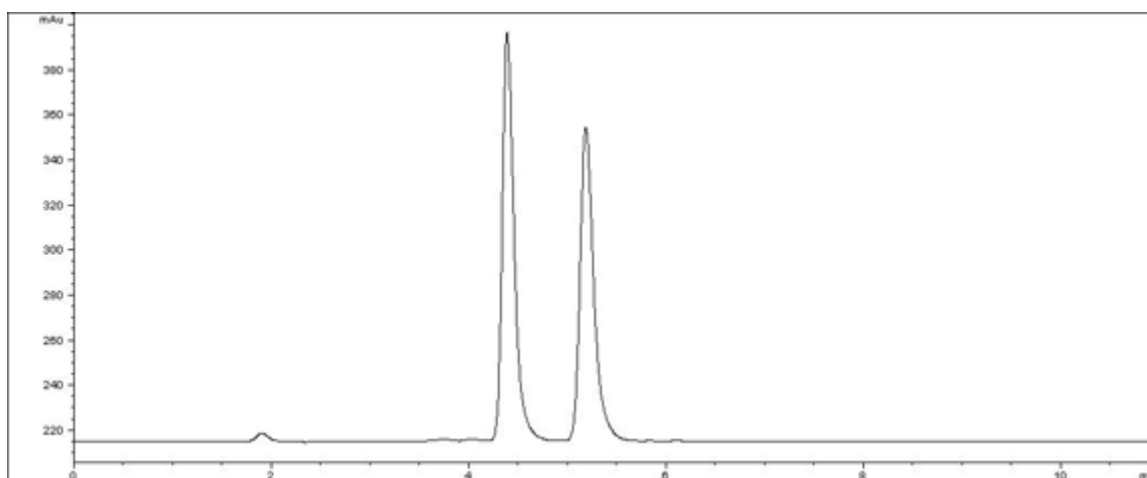

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® ZWIX(+) 150 x 3mm / 3µm
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 98 / 2
Chromatographic Mode	Polar Phase

Flow rate (ml/min)	0.5	Injection amount (mg)	0.0015
Temperature (°C)	25	Pressure (MPa)	10.6
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	4.39	1.43	1 - 2	3.59
2	5.19	1.87		