


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Althiazide

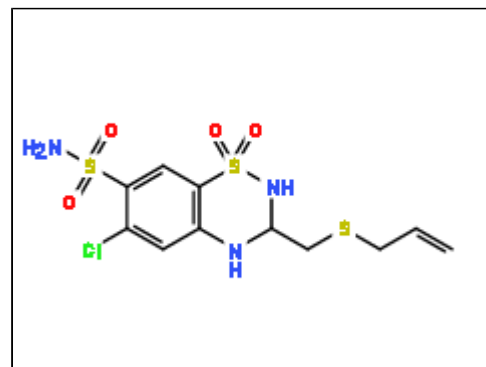

CHROMATOGRAM

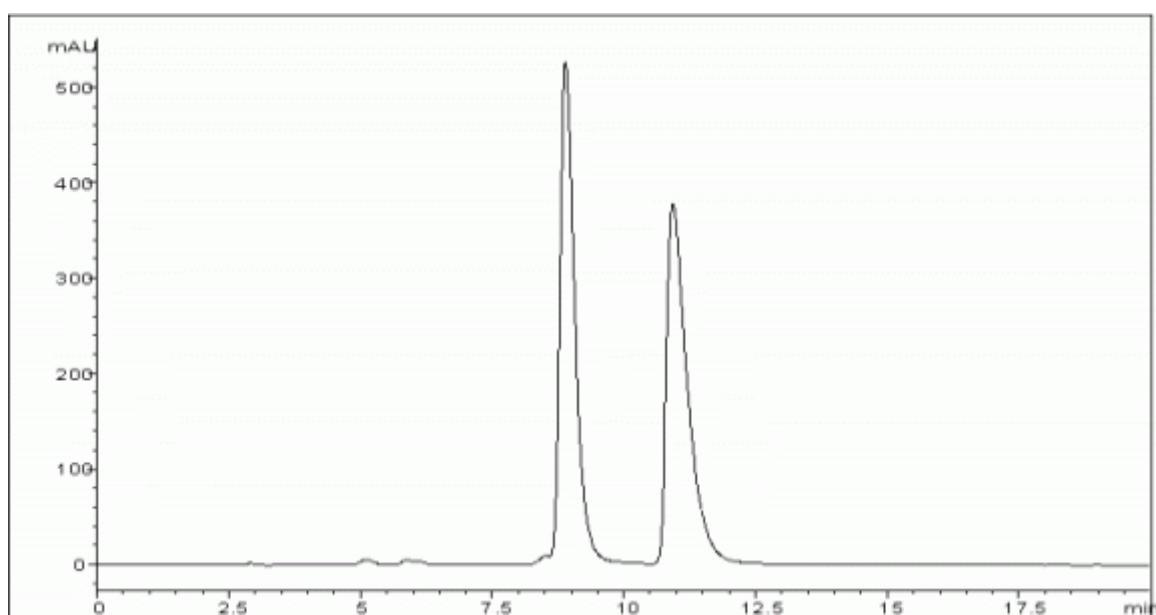

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 75 / 25 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	8.89	1.92	1 - 2	3.44
2	10.93	2.59		