



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Cisatracurium (racemic)
---------------	-------------------------

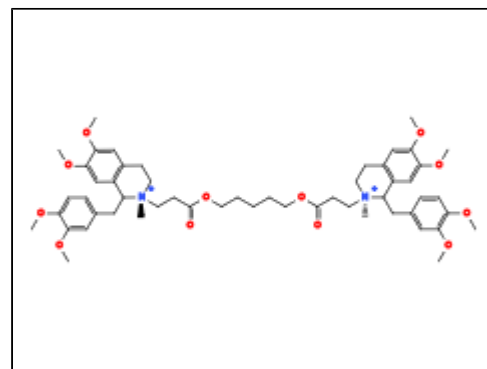

CHROMATOGRAM

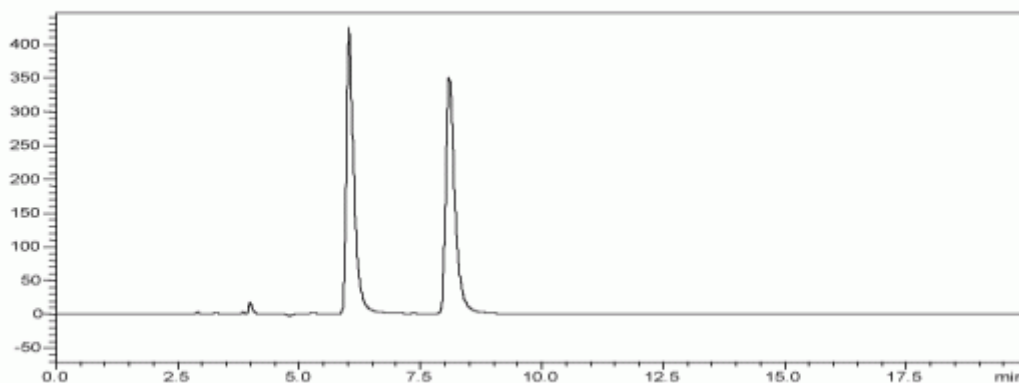

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol / methanol / diethylamine = 80 / 10 / 10 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 240 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.01	1	1 - 2	7.04
2	8.11	1.7		