


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Citalopram

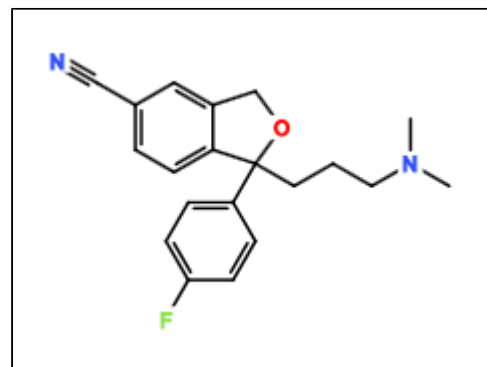

CHROMATOGRAM

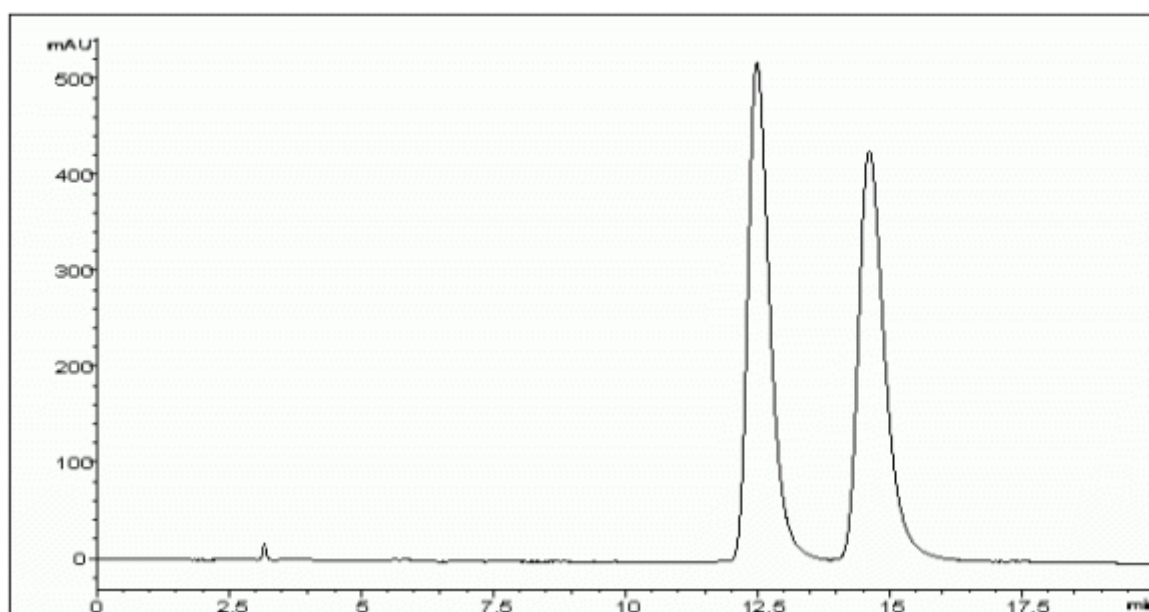

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	n-heptane / methyl tert-butyl ether / diethylamine = 50 / 50 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	12.49	3.16	1 - 2	2.63
2	14.6	3.87		