



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Darifenacin
---------------	-------------

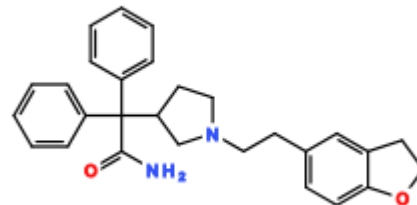

CHROMATOGRAM

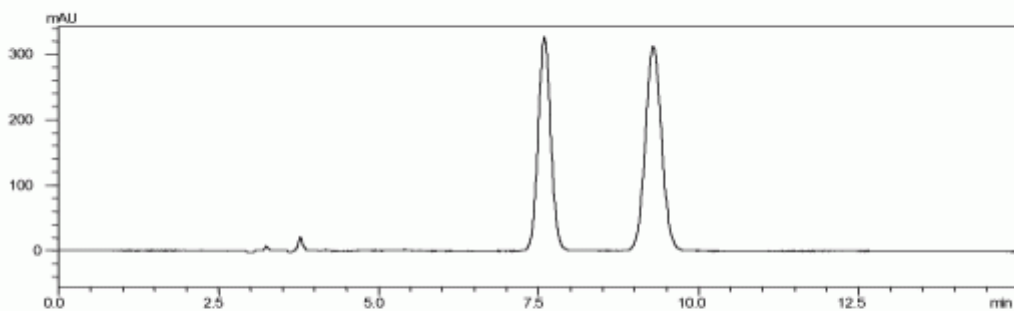

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 50 / 50 / 0.3
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 225 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.61	1.54	1 - 2	1.38
2	9.34	2.11		4