


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Glutethimide

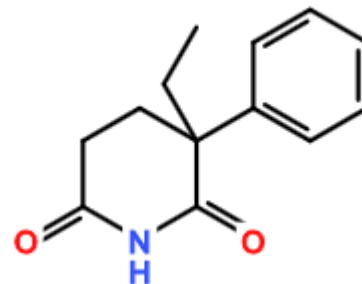

CHROMATOGRAM

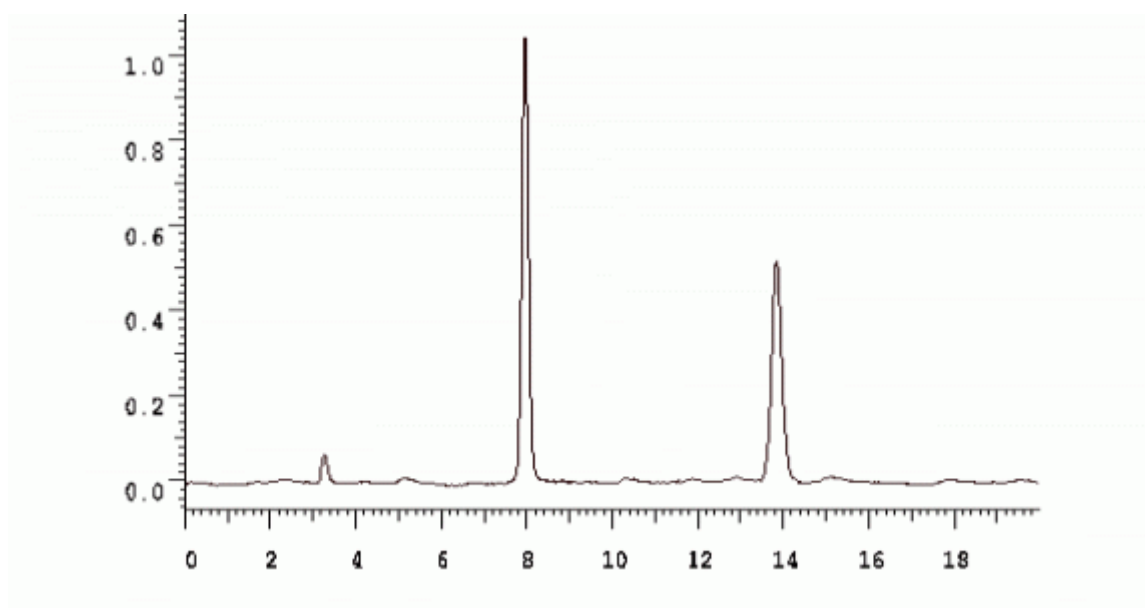

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / dioxane = 70 / 30
Chromatographic Mode	Normal Phase Secondary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.96	1.65	1 - 2	16.31
2	13.85	3.62		