



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Methaqualone
---------------	--------------

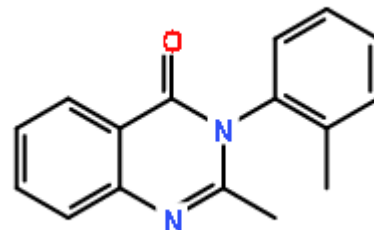

## CHROMATOGRAM

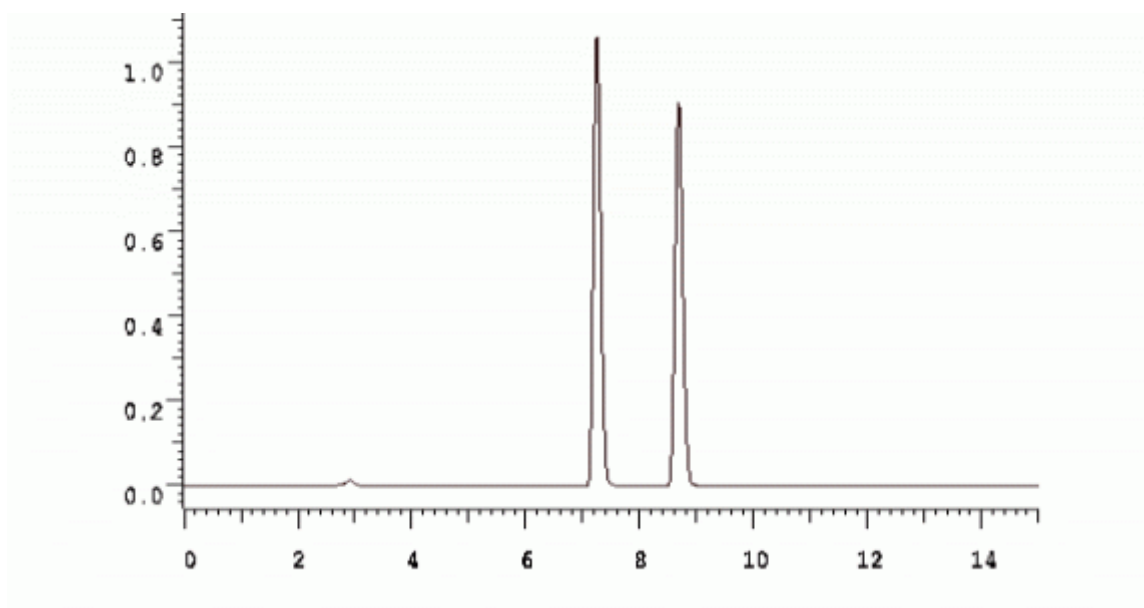

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / acetone = 85 / 15
Chromatographic Mode	Normal Phase Secondary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.01
Temperature (°C)	25	Pressure (MPa)	
Detection	ELSD	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.25	1.42	1 - 2 1.33	5.82
2	8.68	1.89		