



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	N-CBZ-DL-Serine
---------------	-----------------

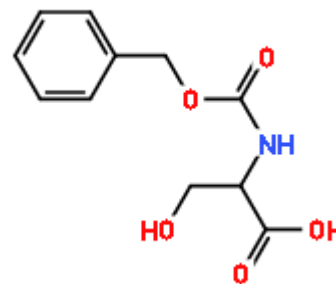

CHROMATOGRAM

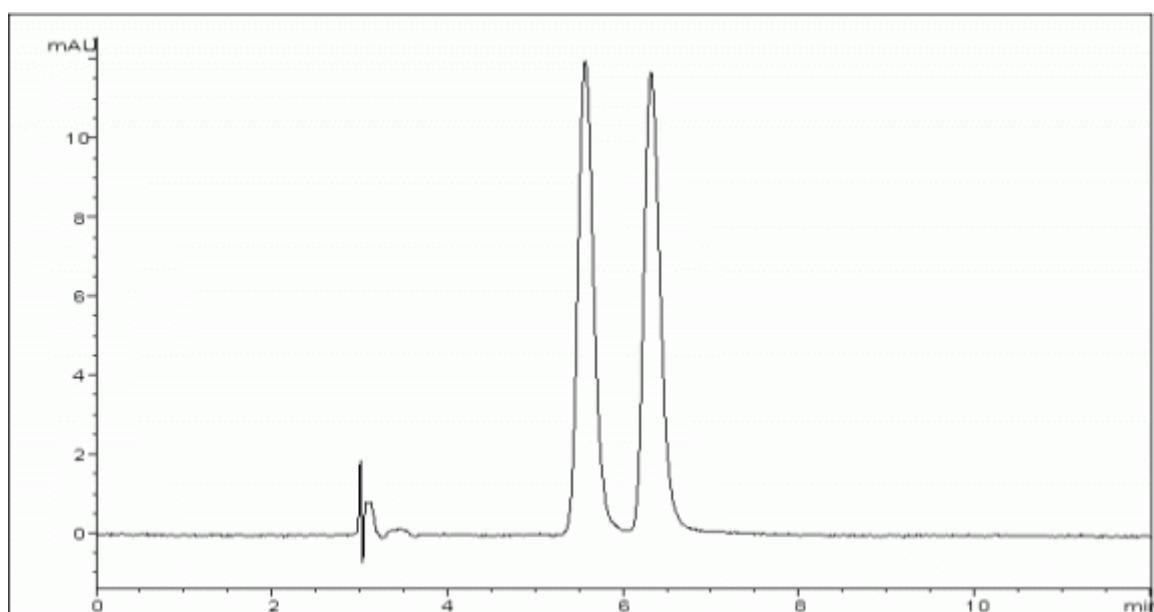

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethyl acetate / trifluoroacetic acid = 50 / 50 / 0.1
Chromatographic Mode	Normal Phase Secondary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.02
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 270 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.57	0.82	1 - 2	2.27
2	6.32	1.07		