


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ornidazole

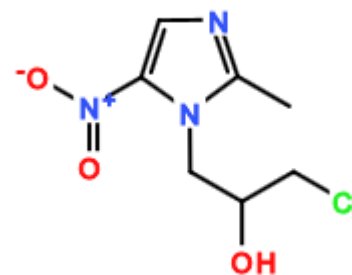

## CHROMATOGRAM

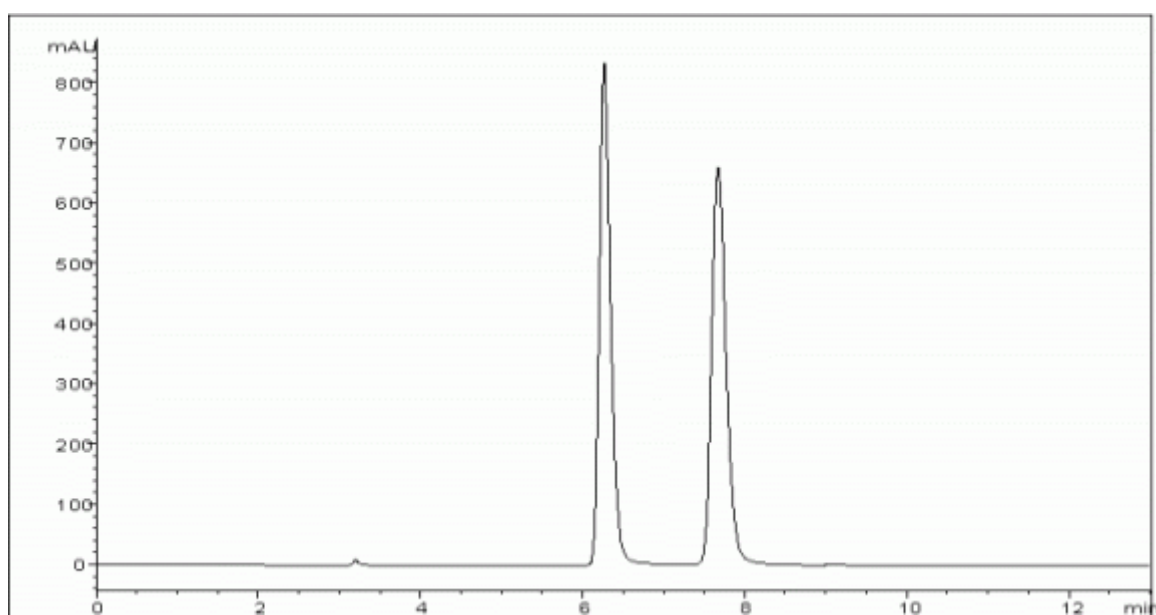

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	methyl tert-butyl ether / ethanolamine = 100 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.02
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.27	1.05	1 - 2	4.96
2	7.67	1.51		