



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Phenylsuccinic acid
---------------	---------------------

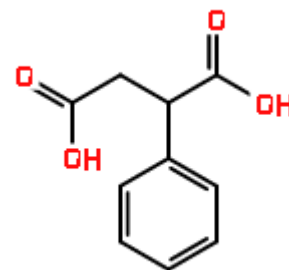

CHROMATOGRAM

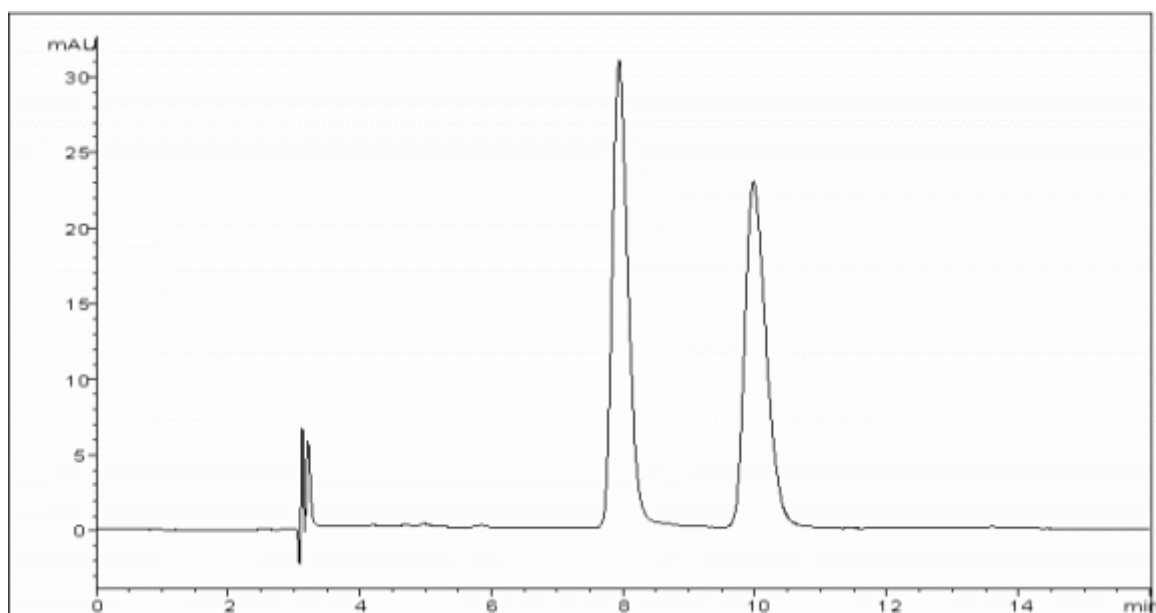

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / methyl tert-butyl ether / trifluoroacetic acid = 20 / 80 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.02
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	7.39	1.59	1 - 2	4.06
2	9.98	2.26		