



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Rosuvastatin (racemic)
---------------	------------------------

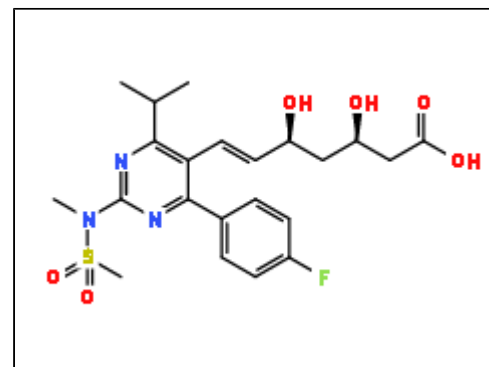

CHROMATOGRAM

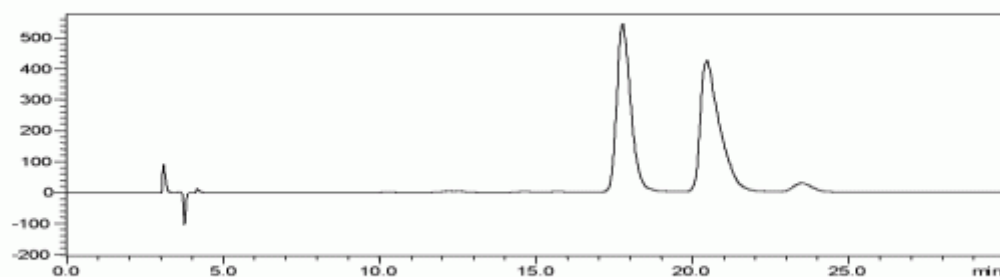

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IB 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / ethanol / 2-propanol / trifluoroacetic acid = 90 / 5 / 5 / 0.3
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 250 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	17.72	4.91	1 - 2	1.18
2	20.44	5.81		2.3