



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Sulconazole
---------------	-------------

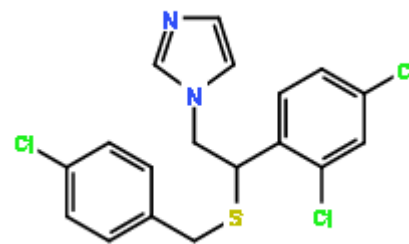

CHROMATOGRAM

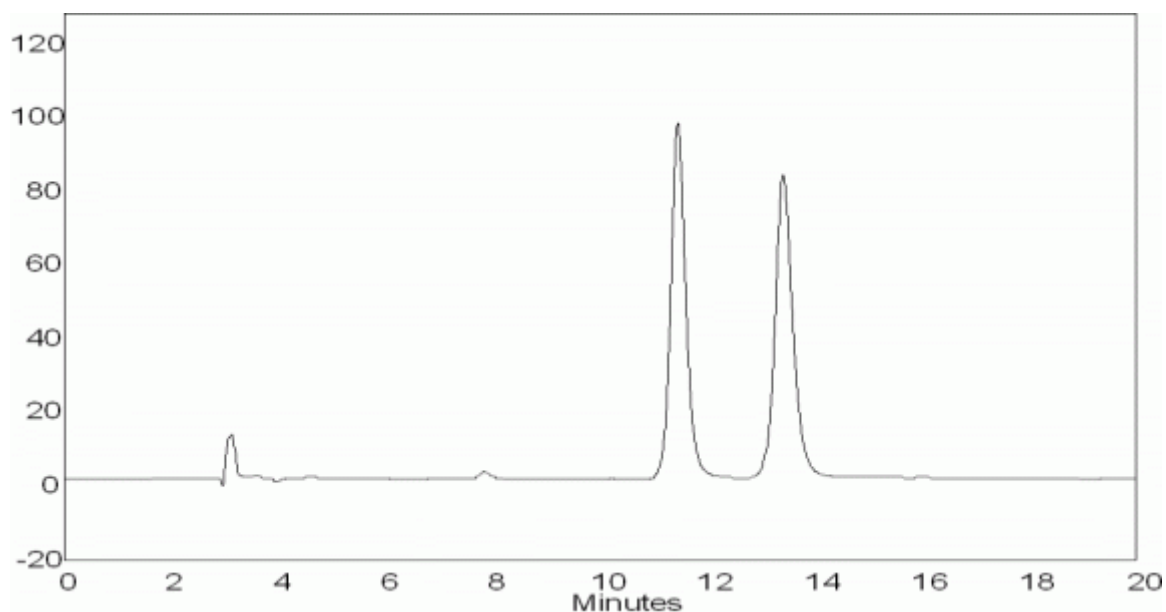

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA 250 x 4.6mm / 5µm
Mobile Phase Composition	methyl tert-butyl ether / ethanol / diethylamine = 95 / 5 / 0.1
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 254 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	11.45	2.66	1 - 2	1.24
2	13.41	3.29		3.55