



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Temazepam
---------------	-----------

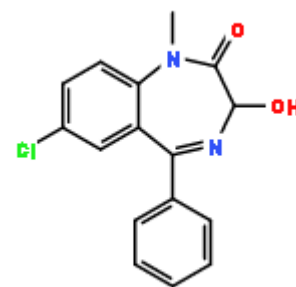

CHROMATOGRAM

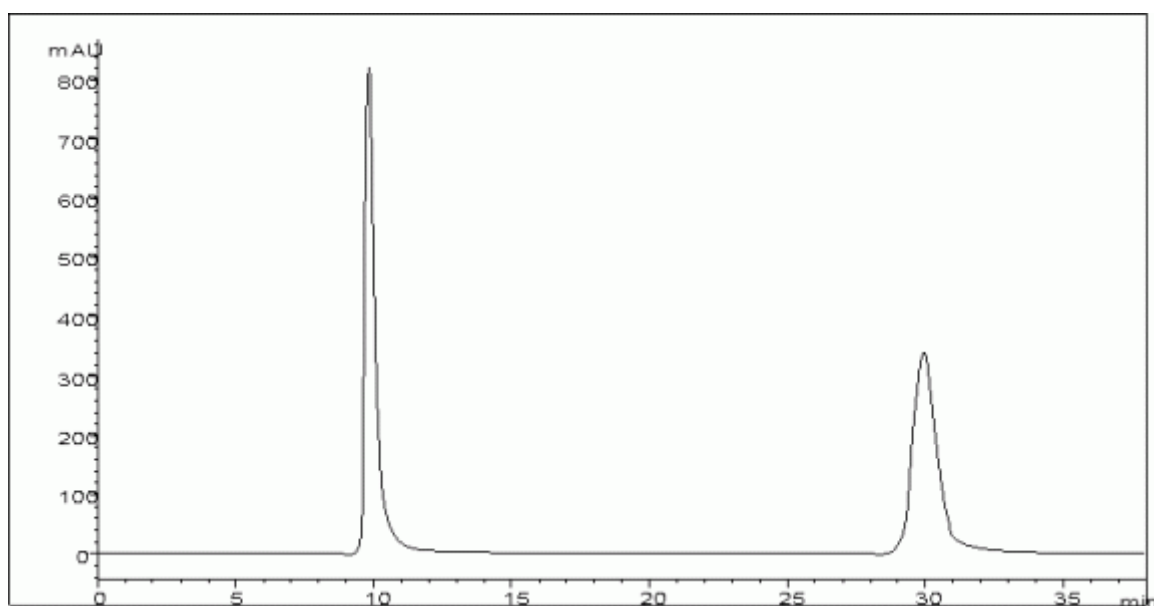

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / 2-propanol = 55 / 45
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.004
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	9.79	2.33	1 - 2	3.94
2	29.91	9.17		19.76