



## APPLICATION RESULTS

## COMPOUND IDENTIFICATION

Compound Name	Tirofiban (racemic)
---------------	---------------------

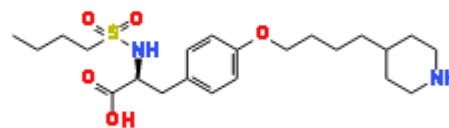

## CHROMATOGRAM

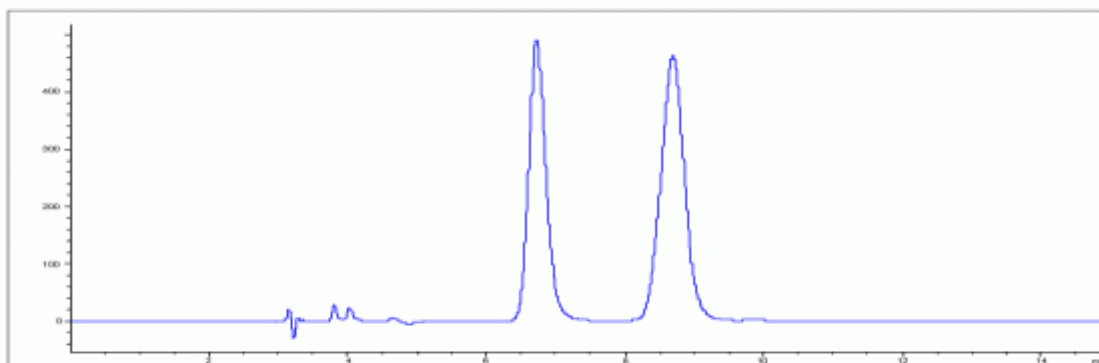

## ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IC 250 x 4.6mm / 5µm
Mobile Phase Composition	n-hexane / methanol / 2-propanol / trifluoroacetic acid / ethylenediamine = 30 / 30 / 40 / 0.3 / 0.3
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	6.72	1.24	1 - 2	1.54
2	8.74	1.91		3.51