



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	Wieland Miescher Ketone
---------------	-------------------------

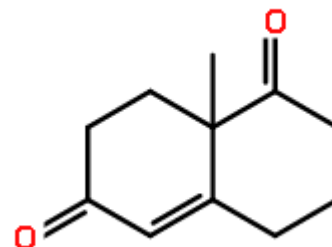

CHROMATOGRAM

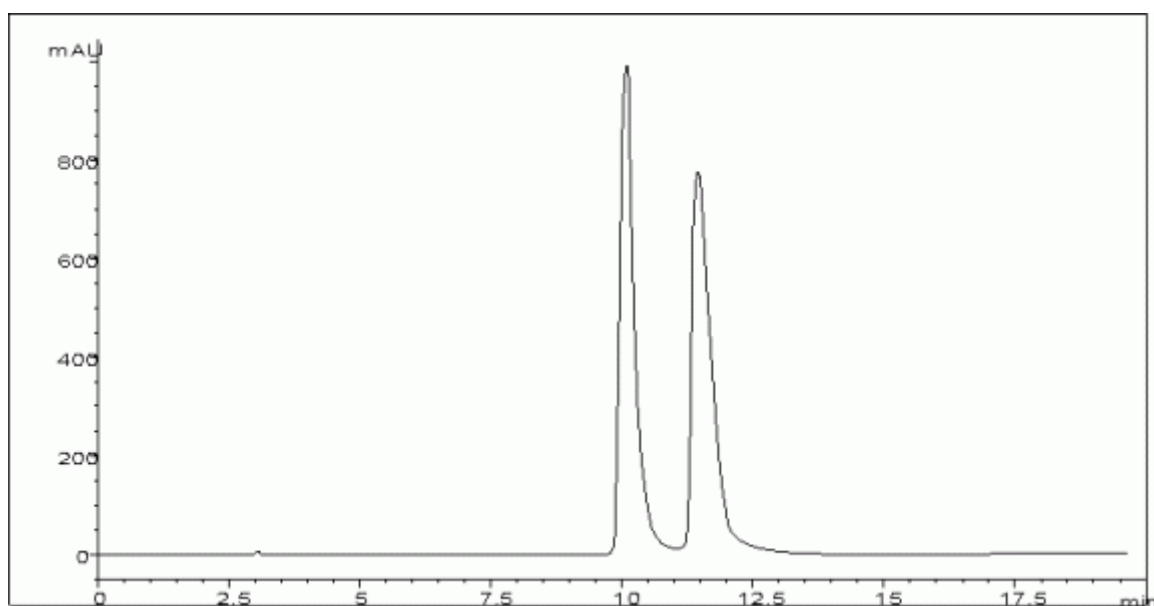

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® IA 250 x 4.6mm / 5µm
Mobile Phase Composition	methyl tert-butyl ether = 100
Chromatographic Mode	Normal Phase Primary screening solvent (Immobilised CSP)

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	10.05	2.35	1 - 2 1.19	2.6
2	11.43	2.81		