



APPLICATION RESULTS

COMPOUND IDENTIFICATION

Compound Name	2-Bromo-1-Indanone
---------------	--------------------

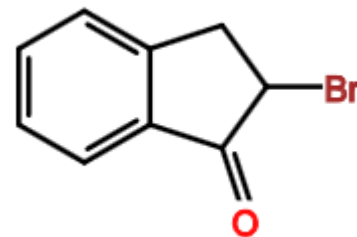

CHROMATOGRAM

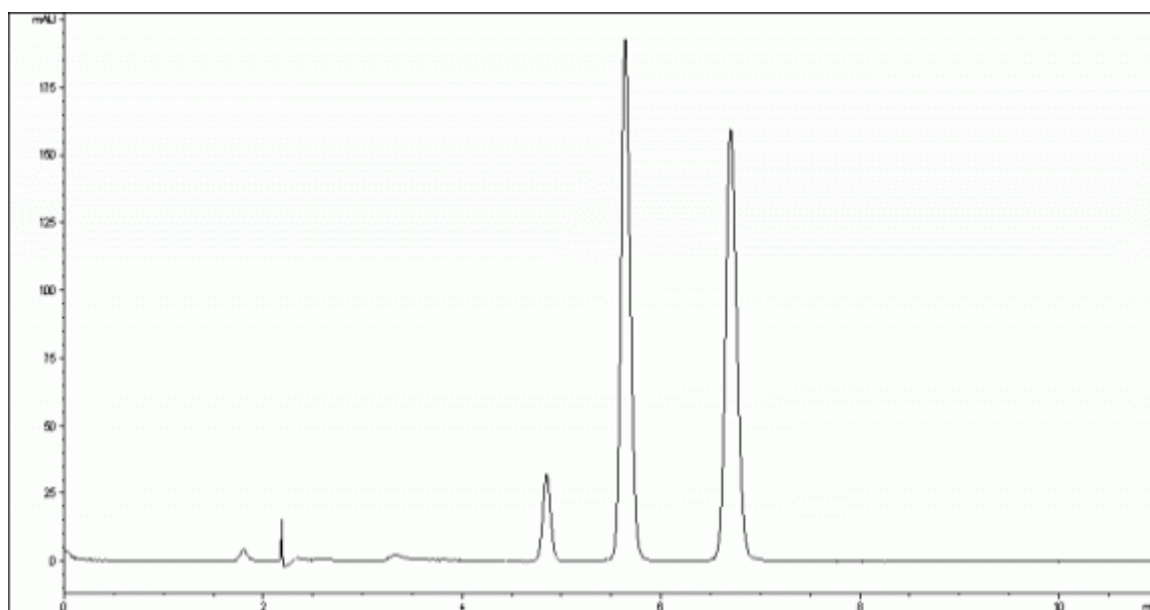

ANALYTICAL REPORT

Column Trade Name Column Size / Particle Size	CHIRALPAK® AY-3 150 x 4.6mm / 3µm
Mobile Phase Composition	n-hexane / 2-propanol = 80 / 20
Chromatographic Mode	Normal Phase Coated CSP

Flow rate (ml/min)	1	Injection amount (mg)	0.005
Temperature (°C)	25	Pressure (MPa)	
Detection	UV 230 nm	Elution order	

No.	Rt	K'	Alpha Value	Resolution
1	5.64	2.19	1 - 2	5.67
2	6.7	2.79		